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Bartosz Machoński

Sonic Electric Toothbrush Toothie X2in1

Instruction Manual

Matters need attention

Forbid

Do not charge in the bathroom. Do not place the charger in a water container above or near a bath, sink, etc.

The battery can not be replaced by the user, only authorized personnel replacement by our company customer service center. Non-professionals are not allowed to replace the product.

Warning

To ensure safety, children aged 3-8, with disabilities, mental or sensory impairment, lack of relevant experience and knowledge under the premise of supervision or guidance use this product. Do not allow children to perform unsupervised activities cleaning and maintenance, the product and its cables should be kept for out of reach of children under 8 years age.

Notice

When using this product for the first time, due to the high speed vibration of the bristles, some people experience tooth irritation or tingling, and even gums may occasionally experience mild bleeding or discomfort, it is normal phenomenon. If this condition still occurs after two weeks of continuous use, please discontinue use of this product and contact your dentist.

If you have recently had oral or gum surgery, suffered teeth problems or two months of extensive oral treatment, please consult your dentist before using the product.



If you feel pain, please stop using this product and contact your dentist.

If the toothbrush is not used for a long time, during charge again, the light does not show, please wait patiently for 6-10 minutes, and then get back to normal.

When charging the toothbrush body feels slightly hot, which is a normal phenomenon. If it feels hot, please turn off the power immediately and bring it to our company authorized customer service center to check.

This product is for personal use and cannot be supplied to dentists and more patients

This product is only for indoor use, do not place the toothbrush directly exposure in the sunlight. It should be stored at room temperature or dry ventilation place

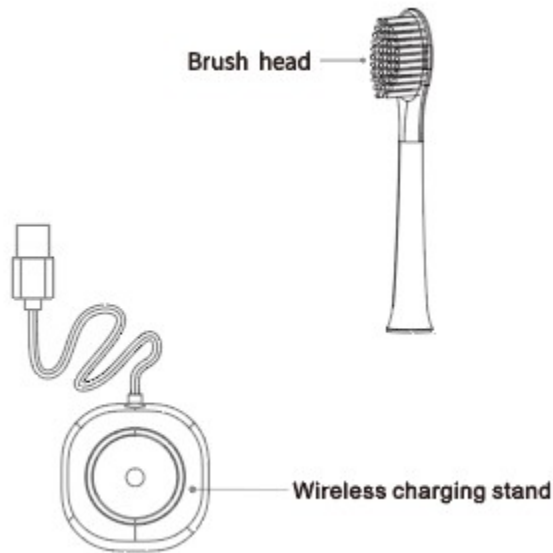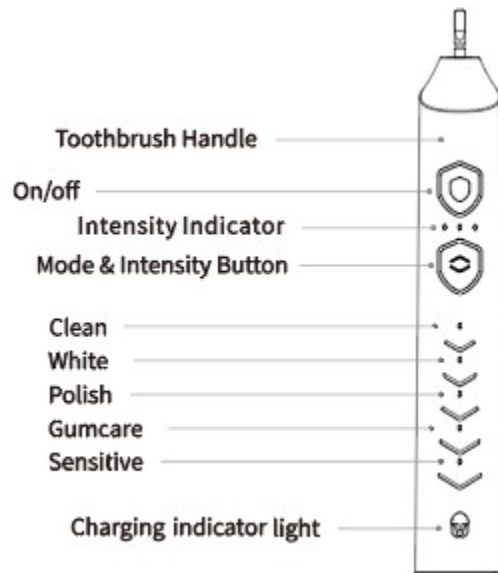
Reminder


Product, accessories, user interface and other illustrations in the manual are schematic diagram. for reference only. Due to product updates and upgrades, products may be slight differences between the real product and the schematic diagram. Please refer to the real products.

Use guide

Step 1: Ensure that the brush head and handle are properly assembled.

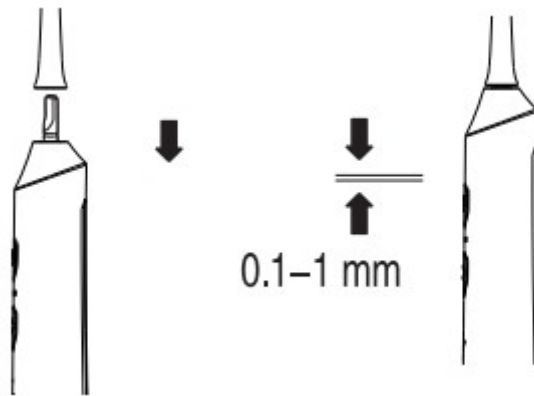Step 2: Wet the brush head and squeeze the toothpaste onto the bristles.

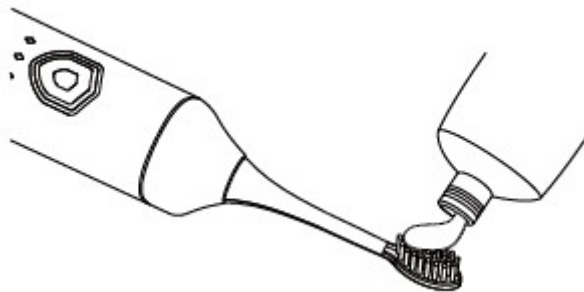Step 3: Power on / mode & intensity switching / Power off.

- Power on: Press the power on / off key;
- Mode switch: When the machine is on, press the mode button and the toothbrush will switch between different brush modes in order;
- Mode strong and weak switching: In the startup state, long press the mode key for 1.5 seconds, the mode work will press weak, medium and strong switching.

Choose your preferred brushing mode:

This toothbrush comes with 5 brushing modes (each mode with three intensity- strong, medium, low;)

Press and hold for 1.5 seconds to switch working intensity:

- Low
- Medium
- Strong

Touch to switch working modes

- Clean: suitable for most people
- White: suitable for smoke stains, tea stains and other black and yellow pigment caused tooth discolouration crowd;
- Polish: suitable for people with dental calculus and tartar
- Gum care: suitable for people with bleeding gums
- Sensitive: suitable for gum and dental sensitive population

Power off

- In working state, lightly press the power button to power off.

- Memory function, after 3 seconds of power on, if the power is turned off, the toothbrush will remember the mode and power level used before the power off, as the default brushing mode after the next power on.
- After each shutdown, the remaining power is displayed by the indicator light of each mode being on for 2 seconds.

Step 4: Brush your teeth correctly

- Place the bristles at a 45 degree Angle from the gum line, starting in a circle move slowly carefully rub the toothbrush against the part of between root and teeth.


- Keep moving the brush head to clean the whole tooth.

- During brushing, Vibration effect of toothbrush will short pauses every 30 seconds, prompt to change the brushing position.
- 2-3 minutes is the best time to brush your teeth, automatically shut down after 2 minutes; if you still need to brush your teeth, you can turn it on and to continue using it.

Step 5: Clean the brush head.

- After brushing your teeth, clean the bristles of the brush head under the facet clean water.
- After cleaning the brush head, please place the electric toothbrush in a dry and cool place.

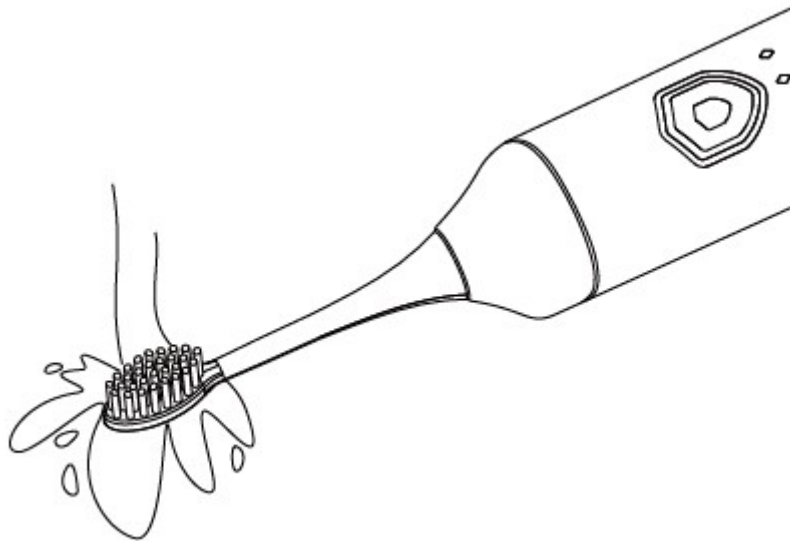
※The dentist suggests changing the brush head every three months.

Maintenance guide

Charging

The red flashing of the toothbrush charging indicator indicates the electricity of the toothbrush. The amount is insufficient, please charge it in time.

Charging method showing as below

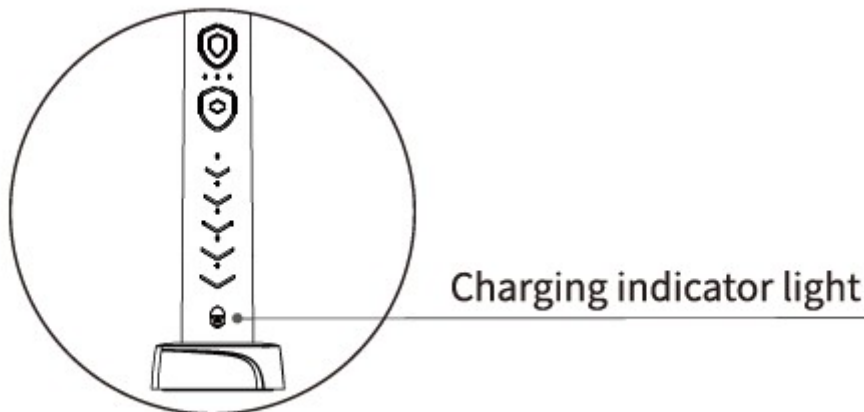- Charging indicator green breathing blinking: charging;
- The charging indicator is green and on the battery is fully charged;
- When the battery is low (the charging indicator flashes red or cannot be turned on), put the sonic electric toothbrush on the wireless charging stand, and connect to the power supply to fully charge for about 12 hours.

Note: The original wireless charging seat is dedicated to charging acoustic electric toothbrush D36, and does not support fast charging for other devices.

Basic data

Wireless charging stand
Rated voltage:5V
Rated power: 2V

The toothbrush
Rated voltage: 3.7V
Rated power: 2W

Charging time: about 12 hours
Water proof grade:IPX7
Origin: China
Execution standard:
GB4706.1-2005.GB4706.59-2008

The charging stand cannot be washed or soaked with water.

Trouble shooting

Problem	Possible cause	Solution
Unable to start the product	Low power	Charging your electric toothbrush
BH cannot insert handle	Foreign material in brush head insertion hole	Remove the foreign material
	Handle driving shaft with foreign material	Clear foreign material
	Unmatched brush head	Purchase official Matching brush head
Short service time	Overuse	Charging to toothbrush
Brush head vibration is less powerful	Brush head and handle assembly improperly	Clear the foreign material reassemble brush head
	Low battery power	Charging to toothbrush

Hazardous substances

Content description

Accessories name	Pb	Hg	Cd	CrVI	PBB	PBDE
Plastic shell	o	o	o	o	o	o
PCB	o	o	o	o	o	o
Charging cable	o	o	o	o	o	o
Packaging material	o	o	o	o	o	o

The environmental protection period of each components: circuit board assembly: 10 years. Power cable: 10 years.

The environmental friendly year identified by this product is no toxic and harmful substances leakage during the period and use the product under the conditions stipulated by this manual.

○ Indicates that the content of this toxic and hazardous substance in all homogeneous materials of certain accessory is within the limitation stipulated by SJ/T11363-2006.

X Represent the content of toxic and harmful substances is out of limitation stipulated by SJ/T11363-2006 at least in one homogeneous material of certain accessory.

 Environment – friend use period

The icon represent 10 years environmental-friendly use period starting from manufacturing date i.e.no toxic and harmful substances leakage and mutation under normal use during the period, safely in use and will not cause harm to environment and human body.

The icon covers all accessories included in sonic toothbrush.

Take the product to an collection point for recycling after 10 years environmental-friendly use period.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

CE declaration

Malis B. Machoński sp. k. Declares that the Maltec device sonic toothbrush meets all the requirements of the EMC 2014/30/EU and meets the following standards:

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