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
Bartosz Machoński

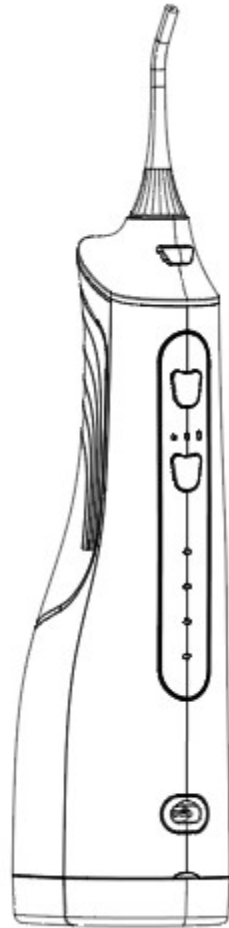**Dental Irrigator
Wireless Toothie X50
Instruction Manual**

CE FC FDA 10 

Statement

To ensure proper and safe use of the product and to avoid damage caused by incorrect operation, please read all instructions and considerations about the product in the manual carefully before using the product.

Thank you for choosing our portable water flosser! Read this manual carefully before using the product. Correctly grasp the use method of it to give full play to the excellent performance of the product, achieving better use results.

Please keep this manual properly and follow it closely so that we can give you better service in the future.

Product Model: WT-288

Operating Mode: Normal/Pulse/Soft/Snipe

Input Voltage: 5V

Input Power: 3W

Full Charge Time: 6H


Waterproof Class: IPX7


Water Tank Capacity: 320ml

Battery Specification: 18650 lithium battery 3.7V, 2600mAh

Execution Standard: GB 4706.1-2005, GB 4706.59-2008

Product Application

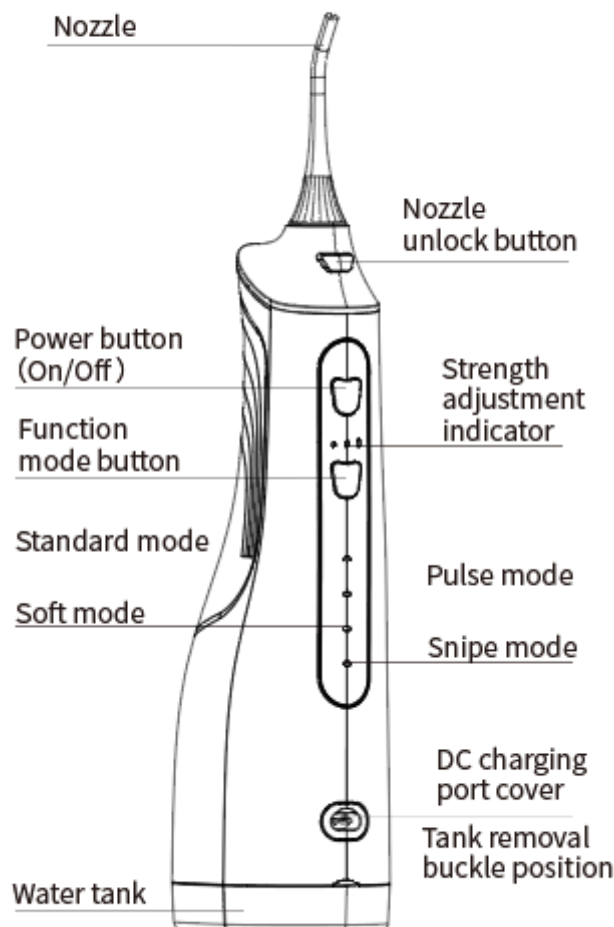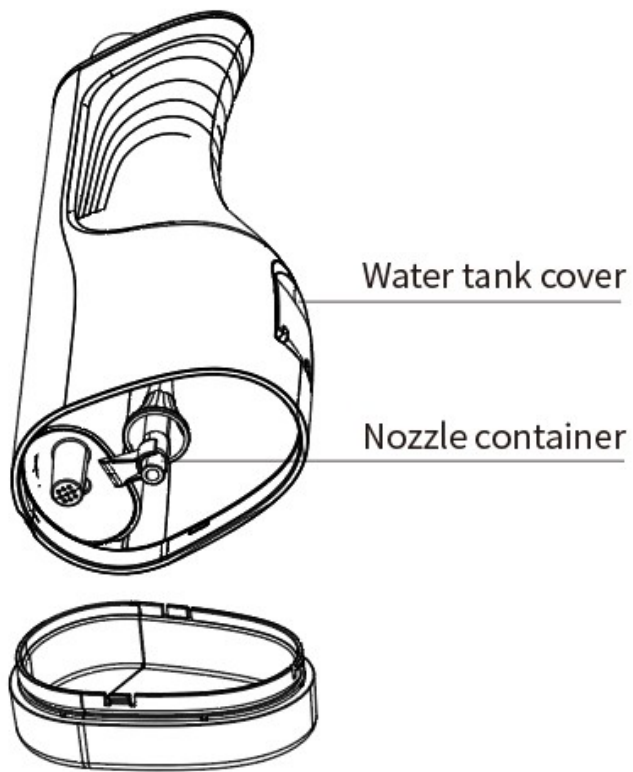Product manual

DC charging cable

Standard nozzle

Orthodontic nozzle

Periodontal pocket nozzle

Illustrations of products, accessories, user interface, etc. in the manual are schematic diagrams for reference only. Due to product updates and upgrades, there may be slight differences between the actual product and illustrated one, please be subject to the actual product.

Safety Precautions

- For users with wounds in the mouth or tongue, periodontal disease or those who have undergone oral surgery in the last six months. Please use the product under doctor's guidance.
- Children under 8 years old shall avoid using the product. When children over 8 years old and people with special needs use the product, using, cleaning and maintenance of the product must be under the supervision and guidance of their guardian.
- Please place this product in a place where children and babies can not touch. This product can not be used as a toy.
- Do not use wet hands to plug and unplug power adapter and charging cable, otherwise it will lead to electric shock.
- Do not drop or insert foreign objects into the charging connector, otherwise it may cause electric shock or short circuit.
- Do not damage, modify, forcibly bend, pull or twist the charging cable, and the USB cable should be kept away from overheated surfaces.
- Before charging, make sure charging connector and charging cable are dry. Forbid charging in humid environment; otherwise it will lead to short circuit.
- Do not use the product for purposes other than oral cleaning.
- Please use a charger or power supply with a USB output port that meets Chinese standards, with output voltage of 5V, output current of at least 1A and output power of no more than 5W.
- Do not use the product without installing jet.
- Please do not fill the tank with hot water above 40 C and clean mouth to avoid scalding.
- Do not spray water directly under your tongue, throat, ears, eyes and other parts of your body, for the pressure of water may cause injury to these parts.
- Do not spray water directly on the gums/periodontal pockets as this may cause bleeding of the gums.
- Do not press the jet hard on teeth, gums or orthodontic braces to avoid damaging gums and jet.
- Do not place the product and accessories in liquids over 40 C corrosive liquids, diluent and alcohol to avoid damaging the product.
- Before use, cleaning and maintenance, please unplug power adapter and charging cable, and keep the product in power-off state.
- If it is used with mouthwash or dentist-recommended solutions, be sure to wash the appliance thoroughly with water after use to avoid contaminating water tank or clogging water pipes with residue.
- After use, please use a piece of dry cloth or dry paper towel to dry appliance moisture, and place it in dry place at normal temperature.
- Do not place the product in hot or humid place and avoid direct sunlight.
- The product has a built-in non-replace-able rechargeable battery. Please do not put it into the fire, heating, or charging, using or placing at a place of high temperature; otherwise it will lead to overheating, fire or explosion.

- Do not disassemble or modify the product, as this may result in fire or personal injury.
- Please unplug power adapter and charging cable when not in use for a long time; otherwise it will lead to electric shock or electrical fire due to decreasing insulation performance. If you find any damage to the product including USB cable and jet, please stop using it.
- Do not use non-original or non-recommended jet and accessories with the product.
- Do not use iodine, salt, bleach, or insoluble essential oil in the product. Use of these will reduce the performance of the product and shorten service life of the product.
- Do not place or store the product in a place where it may easily fall into the bathtub or sink.

Product Application

Charging

- Please use DC charging cable of the product to charge it before first use. It takes about 6 hours to be fully charged.
- Open charging cover. Check to ensure that the product, charging connector and hands are in a dry state. Insert the DC charging cable into charging connector, link power adapter that meets rated voltage (need to provide by yourself) and switch on power.
- Please use a charger that meets Chinese standards or power supply with a USB output port, and make sure its output voltage is 5V and output current of at least 1A.

Nozzle

Nozzle installation/removal

- Insert the nozzle into the nozzle mounting hole, and the installation is complete when you hear a "click".
- Long press the nozzle unlock button and pull it out upward to remove.
- After the nozzle is removed, the water tank can be removed, then the nozzle can be put into it to container for being portable.

Notes

Water tank can be removed: do not press the nozzle unlock button during use to avoid accidental injury from nozzle ejection.

Nozzle Selection

- Standard nozzle: effectively remove food debris, tartar and plaque from teeth and grooves to prevent tooth decay and dental calculus.
- Orthodontic nozzle: wipe and clean the hard-to-reach parts such as tooth sockets or braces. Periodontal pocket nozzle: effectively remove plaque from gingival lines and crevices.

Power On/Pause button

Press On/Pause button, indicator lights up, switching to current mode, press On/Pause button again to start motor, and then short press again to pause the motor work. In Snipe mode (mode light and three intensity lights up), press and Snipe, long press and long Snipe, Auto power off after 10 seconds without any operation in the power-on state.

Mode Button

- In power on state: short press mode button, it switches back and forth in Standard-Pulse-Soft-Snipe mode, all three modes run for 2 minutes, and then automatically power off. Button has no effect in power off.
- Long press the mode button for more than 1.5 seconds for mode strength adjustment, from weak, medium, and strong(1s interval switching in adjustment without letting go)
- Intelligent working memory function: after 3 seconds of working state in power-on, the current working mode will be memorized after auto power off, and next time directly enters this mode.

Low battery reminder

When the battery voltage is lower than 3.2V during operation and startup, the LED (3 strength adjustment lights) flashes (2HZ) as a reminder.

Charging instructions:

When the external input is charged, the 4 mode LED lights are displayed in the running mode (starting from the point injection mode light), and the 4 mode LED lights are long and bright when fully charged. The product does not work when the external input is charged.

Reset function:

Within 10 seconds of input charging, press power button 5 times to restore factory settings (all strength lights flash 5 times 1HZ), that is, restore to standard 1, weak strength. Press any button for normal operation.

6 simple steps to enjoy interdental surfing experience

- Install nozzle;
- Open water tank cover, fill the water tank with fresh water and close water tank cover;
- Press mode button to switch and select needed mode;
- Hold flosser upright, insert jet into your mouth, align between teeth gap or the surface of teeth and gently purse your lips.
- Press power button and clean teeth with spraying water.
- After teeth are cleaned, press power button again to stop spraying water and automatically switch off the flosser.

Notes:

- When using the product for the first time, even healthy gums may tingle itch or bleed. This is caused by the gums not being used to stimuli which will usually disappear within 1 to 2 weeks.
- First time users are advised to use weak strength of standard mode for a week before using other modes.

Repair and Maintenance**Daily use and maintenance**

- After daily use, you need to open water tank cover and pour out the remaining water from tank.
- Remove nozzle and store it in a dry and clean place.
- Dry the flosser with a piece of dry cloth.

Daily cleaning and maintenance

- If water tank cleaning is required. Pull water tank out with little force. Separate it from product, and then wash and dry it with water.
- When cleaning surface of the product, please use a damp cloth, clear water or neutral detergent to clean and wipe. Please do not use hot water over 40 C, corrosive liquid, diluent, alcohol to clean the product and nozzle accessories.
- If not used for a long time, please clean the product and leave it to dry before storage.

Advice

- Water tank is cleaned and maintained every 1 to 3 months to ensure the optimum performance of product.
- Replace standard nozzle and periodontal pocket nozzle every 3 to 6 months. If the tip of nozzle is deformed, replace it with a new one as this can lead to malfunction.
- Replace orthodontic nozzle every 3 months. If the tip of nozzle is deformed, replace it with a new one as this can lead to malfunction.

Environmental Protection

When our portable water flosser reaches the end of its life, it may be necessary to remove internal battery in accordance with the laws and regulations of the country or region where it is located. If it is not possible to remove the battery, please hand the product over to a government-regulated recycling centre.

Name and content of hazardous substances in the product

| Name of Components | Hazardous Substances | | | | | |
|---------------------------|----------------------|----|----|--------|------|-------|
| | Pb | Hg | Cd | CR(VI) | PBBs | PBDEs |
| Circuit control component | O | O | O | O | O | O |
| Motor component | O | O | O | O | O | O |
| Plastic component | O | O | O | O | O | O |
| Hardware component | O | O | O | O | O | O |
| Battery | O | O | O | O | O | O |

| | | | | | | |
|-------------|---|---|---|---|---|---|
| Power cable | O | O | O | O | O | O |
|-------------|---|---|---|---|---|---|

This form has been prepared in accordance with the provision of SJ/T11364.

O: means that the content of hazardous substance in all homogeneous materials of the component is below the limit specified in GB/T26572.

X: indicates that the content of hazardous substance in at least one homogeneous material of the component exceeds the limit specified in GB/T26572

Types of Failures and Treatment

| Performance Failures | Possible Causes | Treatment Methods |
|--|--|---|
| Not working | You have just purchased the product or it has been unused for over 3 months. | 1) Charge the product for 6 hours continuously to give it a full charge. 2) Contact customer service. |
| Not working | 1) There is dirt stuck in the nozzle. 2) deformation of the tip of the nozzle. 3) There is dirt in the pump. | 1) Replace the nozzle. 2) Contact customer service. |
| Only working for a few minutes after each charge | 1) Insufficient charging time and low battery; 2) The battery has reached the end of its life cycle. | 1) Fill the water tank. 2) Adjust the angle of use and use it in a more vertical position. 3) Reinstall water pipe. 4) Contact customer service. |

Notes:

To prevent bacterial reproduction, make sure there is no residual water in the unit. If you don't plan to use it for a long time, make sure to dry it as much as you can before storing it. For first use or long unused, the product not to dispense water temporarily may be caused. Our product can be shaken upside down several times after Power-on.

Non-warranty clauses

- Unauthorized repair, misuse, collision, negligence, abuse, fluid ingress, accident, alteration, incorrect use of non-product parts, torn or altered labelling;
- The quality warranty period has been exceeded;
- Damage caused by force majeure such as fire, flood or lightning strike, etc.;
- Damage caused by failure to use or maintain the product in accordance with the product manual;
- Damage to the product and accessories caused by human factors;
- Use of accessories other than our portable flosser accessories, resulting in product failure or accidents.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

CE declaration

Malis B. Machoński sp. k. Declares that the Maltec device Toothie X50 Wireless Dental Irrigator meets all the requirements of the EMC 2014/30/EU and meets the following standards:

EN IEC 55014-1:2021
EN IEC 61000-3-2:2019+A1:2021
EN 61000-3-3:2013+A2:2021
EN IEC 55014-2:2021

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