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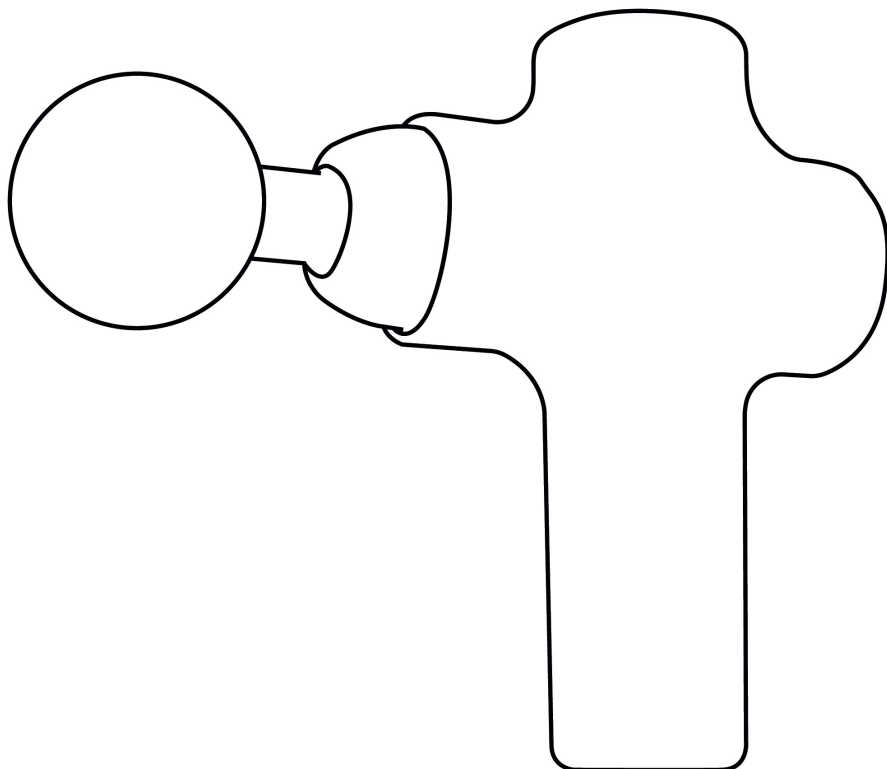
PL

**PIOSTOLET DO MASAŻU
BODYBEAT GO**

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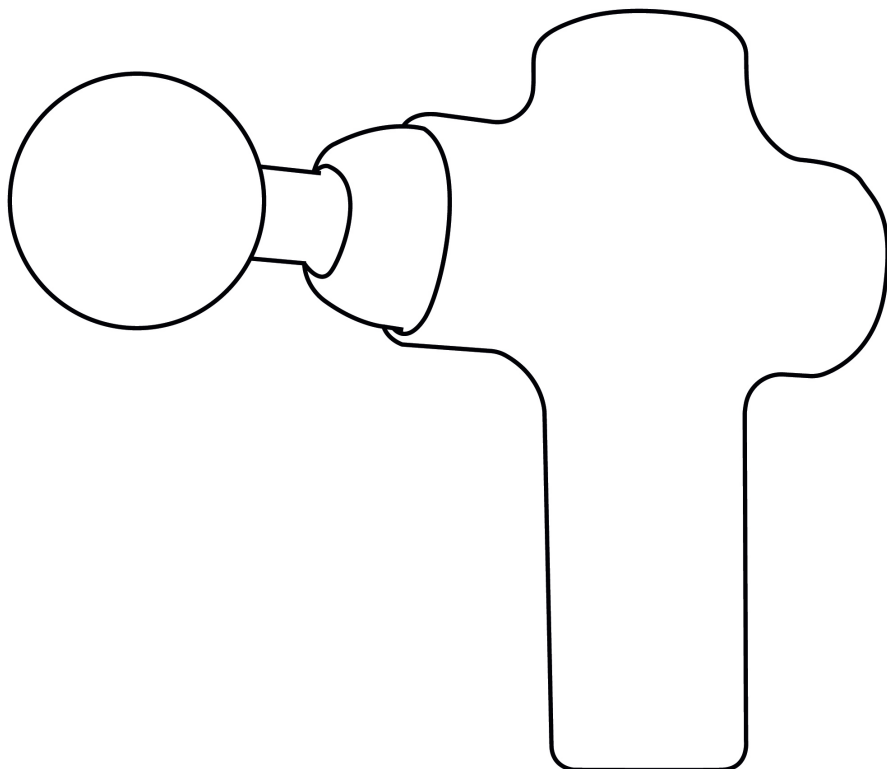
EN

MASSAGE GUN
BODYBEAT GO

**Before using this device, please read this user manual! Follow all safety instructions to avoid damage caused by improper use!
Keep this user manual for future reference.**

SAFETY INSTRUCTIONS

- For use only by healthy adults other than the elderly;
- Only use the product on soft tissues of the body without causing pain or discomfort;
- Do not use on any hard or bony parts of the head or body;
- Whether under moderate pressure or under applied pressure, bruising may occur. Regularly check the massage area and stop using immediately if there are signs of pain or discomfort;
- Can only be used on dry, clean skin surfaces of body parts, gently pressing and moving on the skin;
- Keep fingers, hair, or other parts of the body away from the axis and back of the massage head, as they may be pinched;
- Do not use in areas with less muscle coverage near the bones, such as the tibia, dorsum of the foot, or dorsum of the hand;
- Do not use in or around sensitive areas such as eyes, teeth, breast augmentation, reproductive organs, prostheses, etc;
- Do not use the product in any possible fracture area;
- If the pain is persistent, immediately stop using and consult a doctor. This product is not a medical device and is not recommended for medical use
- Not intended for medical use in hospitals;
- Before use, please check the host and connecting wires. If there is any damage, it must be replaced by the manufacturer or its designated professional personnel to avoid danger;
- Use a massage head that can produce the most ideal effect. Do not place any objects in the
- ventilation holes of the product, do not immerse them in water, and do not let water enter the product
- Please prevent falling or improper use, which may cause permanent damage to the product. Before each use, the operating status of the host, battery, charger, etc. should be carefully checked to ensure that they can meet normal use.
- Do not change or modify the product in any way;
- Please do not let the product run or charge unattended, and do not start the device during the charging process
- Prohibit throwing equipment into fire;

- Do not use the device continuously for more than an hour. Stop the device for 30 minutes before reusing it;

BATTERY

- Do not throw the battery into water or place it in a high temperature and humidity environment!
- Do not charge the battery under fire or extremely hot conditions! Do not use near heat sources such as flames or heaters. Store the battery! If the battery leaks or emits an odor, it should be immediately moved to a safe place and placed separately; Do not remove. Put the battery into the fire or heat it up!
- It is prohibited to short circuit the positive and negative poles of the battery with wires or other metal objects, and it is also prohibited to connect the battery to a necklace, hair clip, or other accessories. It is transported or stored together with metal objects!
- Do not use nails or other sharp objects to pierce the battery casing, and do not hammer or step on the battery!
- Do not hit or throw batteries. Do not disassemble the battery in any way!
- Do not place the battery in a microwave or any pressure container! If the battery emits odor, heat, deformation
- Do not use when discoloration or any other abnormal phenomena occur; If the battery is in use or charging, immediately remove it from the instrument or disconnect the charger and stop using it!

CHARGING

- Before initial use, please charge for approximately 1.5 hours to fully charge the battery;
- Connect the charger interface to the DC charging port at the bottom of the massager, and then connect the charger input to the mobile power supply or mobile phone charger USB output port;
- Please be sure to use a 5V-1A charger or mobile power supply 5V-1A output port;
- In the charging state, the battery indicator light flashes blue, and all four indicator lights turn on blue when fully charged;
- Please ensure that the massager is turned off when charging. When the device is running, charging will automatically protect it, and the battery can be charged at any level,
- Please do not completely drain the battery and charge it in a timely manner;

USAGE METHOD

Preparation before use

- Open the packaging box, check if the host and accessories are intact and complete. If there are any defects, please notify your supplier immediately;
- Ensure that the device battery has sufficient power.

OPERATING STEPS

1. Instrument connection massage head With the device turned off, fully insert the massage head with pressure. When removing the massage head, simply pull it out (note: do not install or remove the massage head while the massager is running, as it may cause mechanical injury).
 2. In the shutdown state, long press the „ „ button on the panel for 2 seconds to turn on the device and enter the first gear. Short press the „ „ button on the panel, the gear will switch between 1-4 gears. When the device is running in 1-4 gears, long press the „ „ button for 2 seconds to stop output, and the instrument will shut down.
- The device automatically stops and shuts down after 10 minutes of operation. If you need to continue using it, you need to restart it
 - LED light indicates the corresponding strike frequency level, with-gear 1 being the lowest and gear 4 being the highest

ABOUT MAINTENANCE

- Clean the surface of the massager with a damp towel, then dry it with a soft cloth. When not using the device for a long time, please store it in the packaging box
- Do not disassemble screws or attempt to disassemble the product;
- To ensure optimal use of the new battery in a short period of time, it is recommended to conduct at least three complete charges and discharges. When the device is not in use for a long time, the battery should be charged once a month

PRODUCT INSTALLATION AND TRANSPORTATION STORAGE STRIP

Product installation conditions

- Relative humidity: $\leq 90\%$;
- Atmospheric pressure: 70kPa-106kPa;

Transport

- Environmental temperature: -20 °C -50 °C;
- Relative humidity: ≤ 90%;
- This product is suitable for transportation methods such as cars, trains, ships, airplanes, etc
- The product should be fixed during transportation to prevent
- Severe vibration and collision
- Rolling and falling, upward, careful, handled with care, and protected from rain;
- The number of layers stacked for massagers is strictly prohibited from exceeding 5 layers for Storage
- The massager should be stored in a well ventilated, non corrosive gas warehouse, 250mm above the ground;


The number of layers for massager stacking is strictly prohibited to exceed 5;

PRECAUTIONS FOR PRODUCT DISPOSAL

Please support the recycling of rechargeable batteries (carp ion batteries).

	The battery used in this product is a specialized rechargeable battery and cannot be used in other products, • Do not charge a rechargeable battery that has been removed • Do not throw it into the fire, as it may overheat excessively • Please do not use nails to pierce, impact by gravity, or undergo decomposition and modification • Do not directly connect the positive and negative electrodes of the battery with metal • Do not store it with necklaces, hair clips, etc • Please do not charge, use, or place it in high-temperature places such as by the stove or in the scorching sun The above situations can all cause heating, fire, rupture, and other phenomena
	Please keep the removed rechargeable battery out of reach of infants and children • If swallowed by mistake, please seek medical attention immediately • When battery fluid leaks, do not touch the liquid with your bare hands and handle it as follows • When battery liquid enters the eye, there is a risk of blindness. Do not rub the eyes directly with your hands, • Immediately wash with clean water and seek medical attention • When battery fluid gets on the body or clothing, it may cause skin inflammation or injury. Please rinse with clean water Wash, if you still feel unwell, seek medical attention immediately • Only products can be decomposed when discarded. Failure to comply may result in fire, electric shock When accidents such as injuries occur and repairs are required (except for battery replacement), please consult the original seller.

MULTIPLE MESSAGE ATTACHMENTS FOR CHOICE

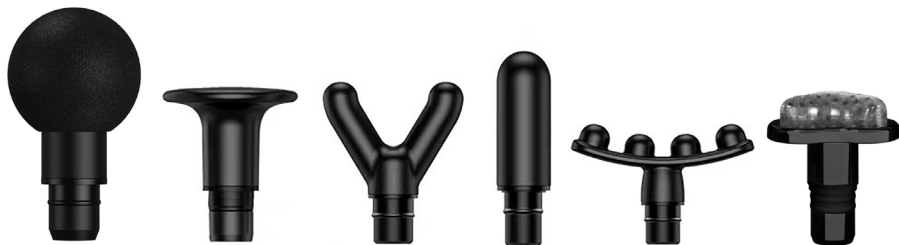TECHNICAL DATA

SYMBOL	113859
INPUT	5V / 2A
OUTPUT	7.4V
SPEED LEVELS	4
BATTERY	Li-Ion 2000mAh
CHARGING TIME	1-2H
USAGE TIME	1,5H
AMPLITUDE OF STROKES	10mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

Li-ion

DECLARATION OF CONFORMITY

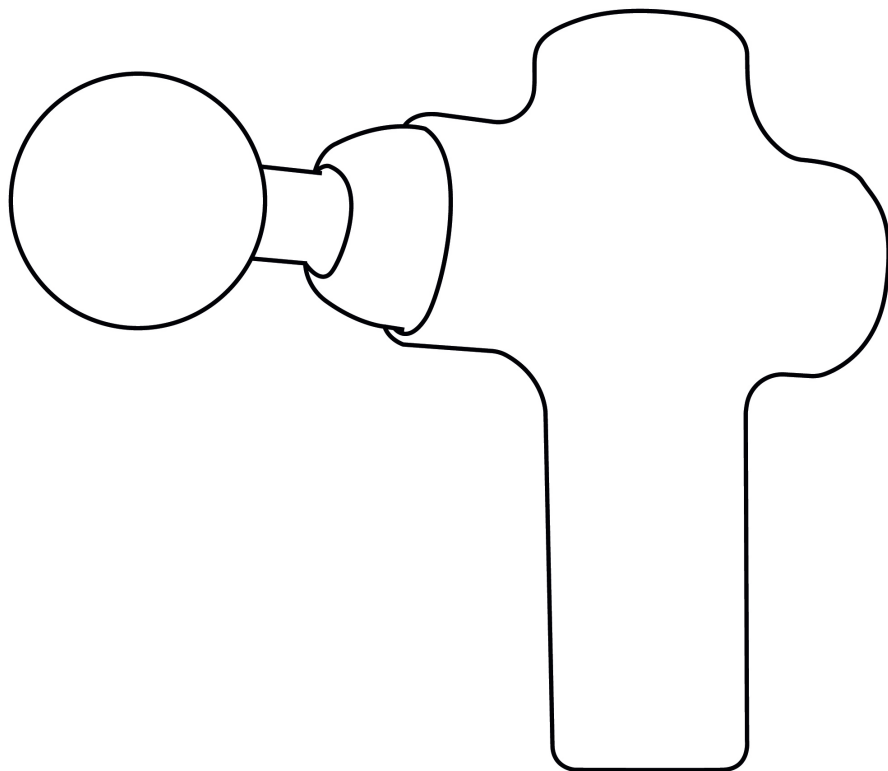
Malis B. Machoński Sp. K. declares that described in this manual Maltec MESSAGE GUN BODY BEAT GO meets all the requirements of the:

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