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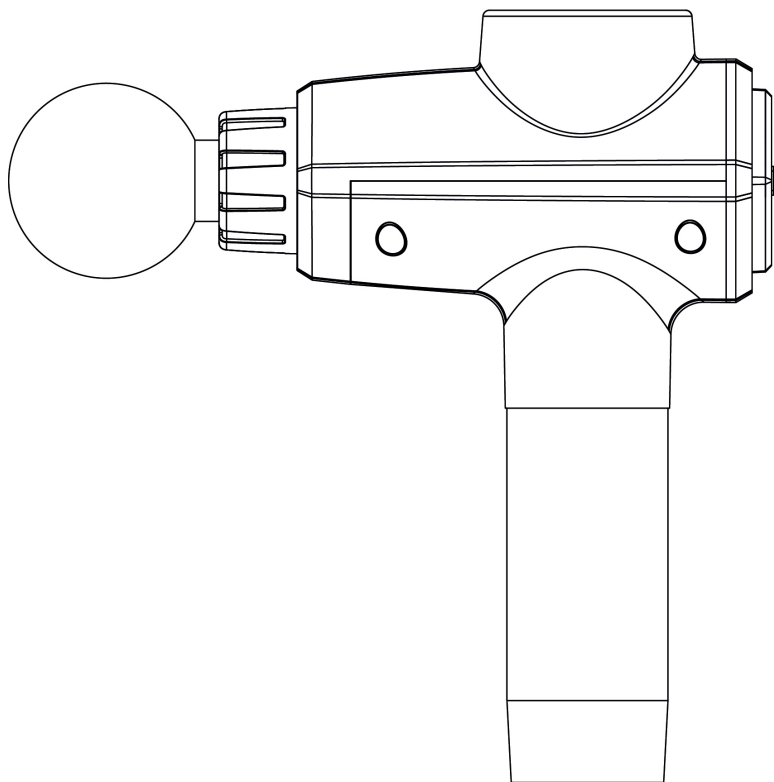
PL

PISTOLET DO MASAŻU
BODYBEAT PRIME

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
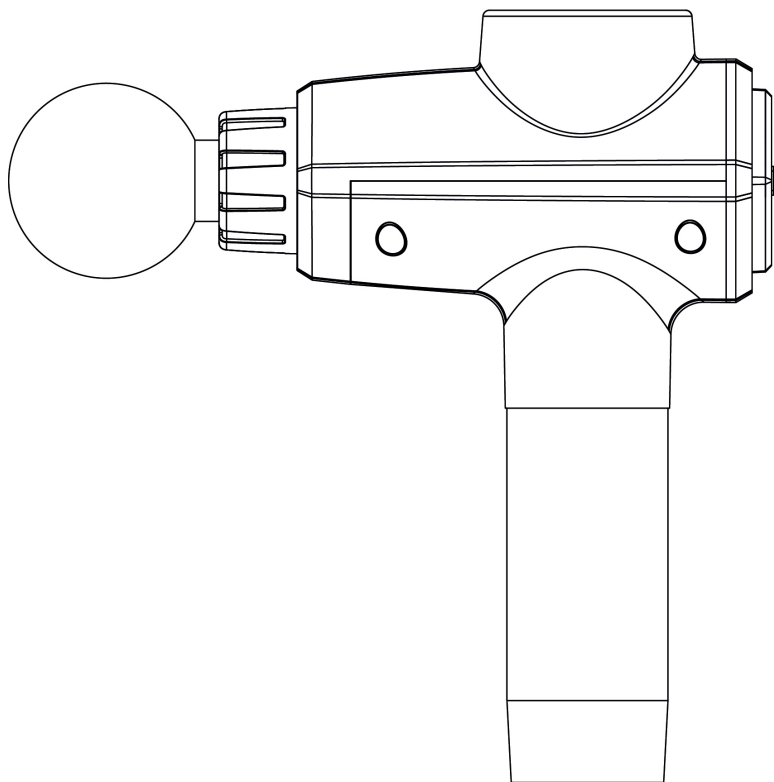
EN

MASSAGE GUN
BODYBEAT PRIME

**Before using this device, please read this user manual! Follow all safety instructions to avoid damage caused by improper use!
Keep this user manual for future reference.**


**1 - round head | 2 - Air cushion head | 3 - spine head | 4 - bullet head
| 5 - flat head | 6 - hand head | 7 - Fan-shape head | 8 -Shovel head |
9 - fork head | 10 - cooling head**

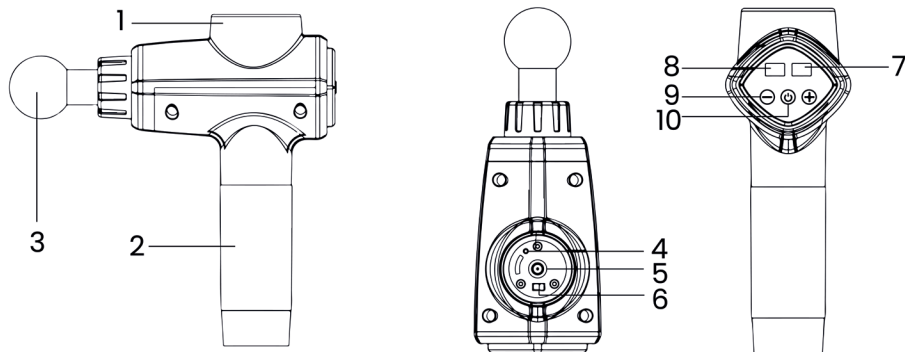1. Brushless motor
2. Rechargeable battery
3. Six massage heads
4. Charging LED
5. Type-C charging jack
6. Power switch on/off
7. LED display for battery left


8. LED display for gears (20 Speed Version)
9. Touch button for speed up/down (20 Speed Version)
10. Touch button for stop (20 Speed Version)

HOW TO USE

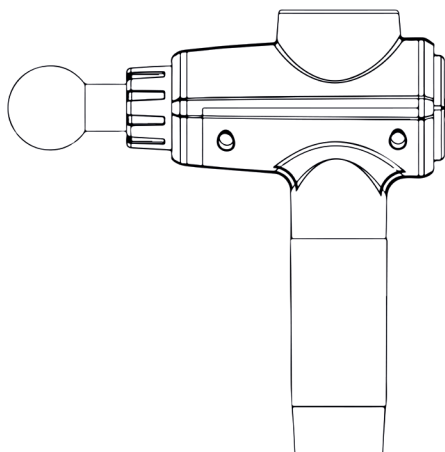
STEP 1

Select a suitable massage head and install

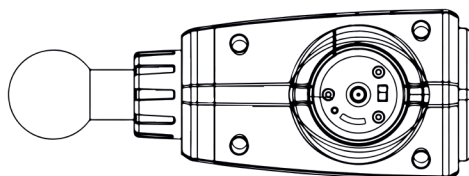STEP 2

Switch on

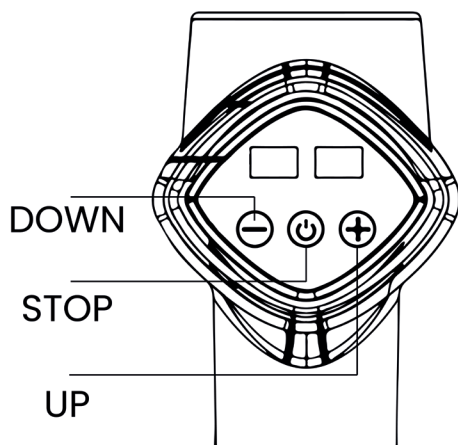STEP 3

(20 Gear Version):

Up/Down: one time touch only changes one gear. Continuous touch to reach lower or higher gears quickly.
Stop: Motor reduces to 0 speed and enters standby status.

NOTES

When you don't use the device and not switch off, 2 minutes, LED display will switch off to save power.

Friendly definition in software:

a. For touch version, press and hold the touch button, can fast change speed to maximum or minimum.

b. For touch version, one button operate can reach any speed, the speed is loop out.

c. In standby status, touch any button reach gear 1.

d. For button version, press and hold the button, can switch to 0 gear.

When you charge the device, when the battery less than 10%, you can not use it

When the battery left only 5%, the device will auto stop to work to protect the battery.

When you charge the battery, the charge LED will flash, when it is full, LED will stop flashing.

And when you charge the battery, battery LED will be on, show the capacity of battery.

When you use the speed you like, next time when you use it will stay in the same speed.

BATTERY INSTRUCTIONS

When the battery is empty it needs to be charged for 2.5–3 hours to be fully charged. The battery life is determined by the gear adjustment. When the battery is charging, the battery power-on indicator light (5) will light up, indicating normal charging. For your safety, the motor button cannot be activated when the battery is charging. Starting the motor during charging is prohibited.

SAFETY INSTRUCTIONS

To reduce the risk of electric shock, fire and personal injury, this product must be used in accordance with the following instructions:

1. Only to be used by adults.
2. This equipment should not directly massage the skin surface to avoid damage to the skin. The massage part should be covered by clean, dry clothes or a clean, dry towel. It should be gently pressed and moved. Massage time on the same part should not exceed 60 seconds.
3. High speed and pressure should only be used on soft tissue parts

of the body, if there is no pain or discomfort. High speed and high pressure must not be used on the head or any bones of the body.

4. Please choose the appropriate massage head, according to the massage needs of the different parts.
5. Frequent high-speed high-pressure use exerts pressure on the same body part. If abrasions, painfulness or incapacity occur, please stop using immediately.
6. Fingers, hair or any part of the body should be kept away from the telescopic rod of the device and the cooling vents of the back motor to prevent entanglement.
7. Do not put any objects into the cooling vents of the motor.
8. Do not immerse the device in water. Do not allow water to enter the device or the motor cooling vents.
9. Do not misuse the device. Dispose of the device correctly after the end of its useful life.
10. Only use the charger provided by the original manufacturer to charge.
11. Please check the equipment and battery carefully before each use.
12. Do not make any modifications to the equipment.
13. Do not run or charge the device unattended.
14. This device may not be disposed of with the normal household waste at the end of its life, but must be taken to a special collection point for the reuse of electrical and electronic devices. The materials used in the device can be recycled. By recycling used household appliances you make an important contribution to the protection of our environment. Ask your local government for the collection points.
15. The device is only intended for massaging the human body. For damage caused by the use of the device for other purposes, incorrect operation or repairs carried out by an amateur, no liability will be accepted.

MAINTENANCE

1. Wipe the surface of the device with a slightly wet towel and dry it with a soft piece of cloth.
2. When storing or traveling for a long time, separate the device from the battery and place it in the device's original case.

WARNING AND CAUTIONS

1. Do not put the device inside water, keep the machine away from

liquids and heat sources, ensure the motor cooling holes are free from dust and debris.

2. Do not remove the screws by yourself.
3. Turn off the device after one hour of operation and continue using it after resting it for 30 minutes.
4. Remove the charger from the battery after charging, and do not use the device during charging.
5. The use of this device for the injured part of the body is prohibited without consulting the doctor.

Please read the instruction manual carefully before using this equipment, use this equipment correctly in accordance with the instructions, our company will not be liable for any damage or loss to the customer caused by the incorrect use of this equipment as explained.

DO NOT USE THIS DEVICE WITHOUT THE APPROVAL OF YOUR DOCTOR IF:

- pregnancy
- diabetes
- neuropathy or retinal damage
- pacemaker wear
- had surgery recently
- migraines
- spinal inflammation
- IUD's
- have had/have health problems

You must consult your doctor first and use it with your doctor's permission.

We are studying the effects of the device on specific medical conditions, this may cause the above taboo symptoms to be not detailed enough, if you are unwell please stop using, follow the advice of a physician or physiotherapist.

We provide warranty for battery, motor, controller of one year, other wearing parts warranty of three months.

EXCEPT WHEN:

1. Damage to this equipment caused by man-made incorrect operation or transport.

2. Unauthorized removal and repair of this equipment.
3. The instruction manual isn't being respected.
4. Overloaded use of this product design parameter.
5. Damage due to the customer's unusually poor storage or maintenance environment.
6. Failure to provide proof of purchase date.

TECHNICAL DATA

SYMBOL	113858
INPUT	5V / 2A
OUTPUT	11.1V
SPEEDS LEVELS	20
BATTERY	Li-Ion 3000mAh
CHARGING TIME	3h
USAGE TIME	4h
AMPLITUDE	14 mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

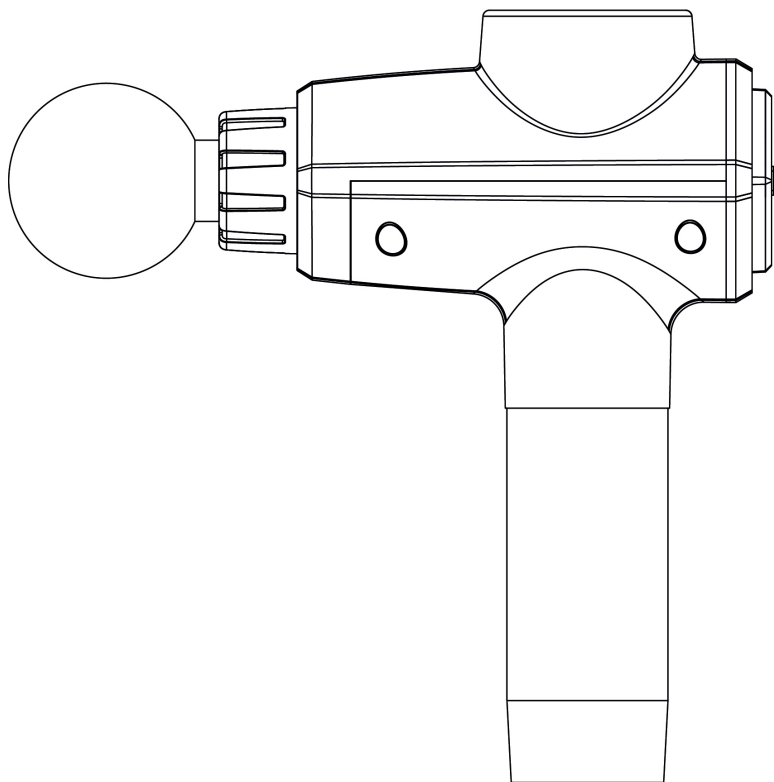
Malis B. Machoński Sp. K. declares that described in this manual Maltec MESSAGE GUN BODY BEAT PRIME meets all the requirements of the:

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