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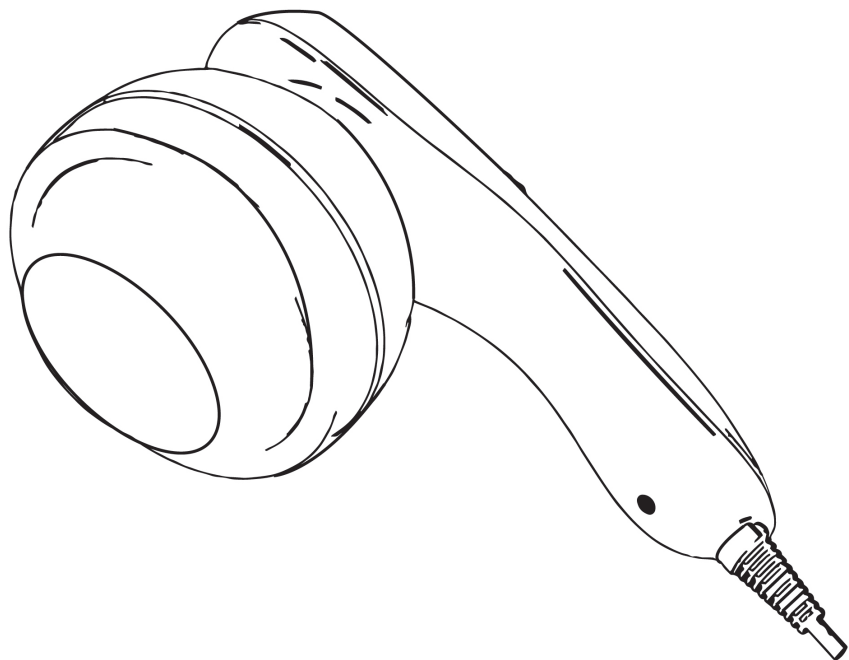
PL

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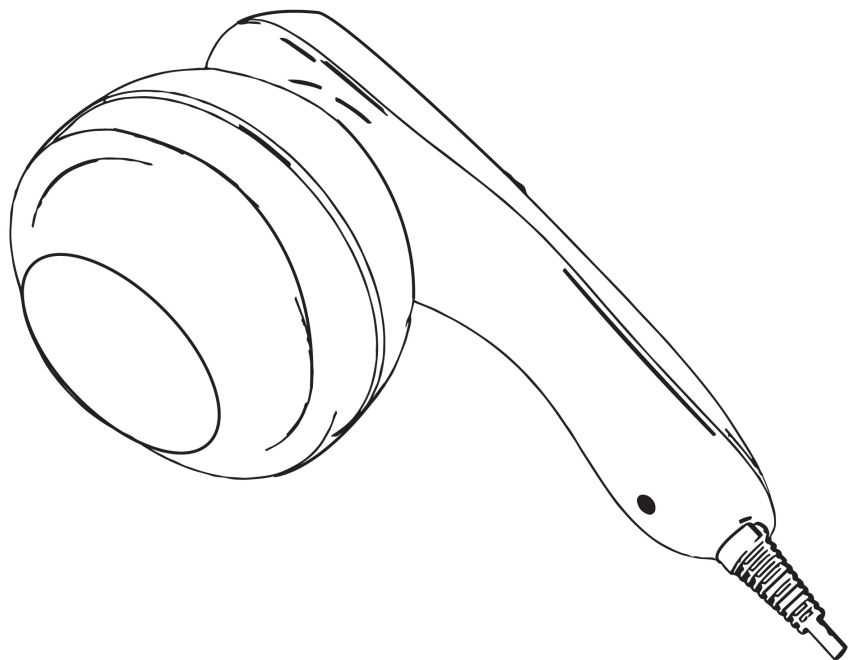
EN

BODY MASSAGER
SKINPULSE

**Before using this device, please read this user manual! Follow all safety instructions to avoid damage caused by improper use!
Keep this user manual for future reference.**

SAFETY INSTRUCTIONS

- Do not use the device if you are pregnant or suffering from heart problems or any physical injuries, and consult your physician if you are unsure.
- Only use the device one hour after meals.
- Keep away from children.
- Do not use the device on the same area of body for more than 15 minutes at a time.
- Do not use the device in an environment that exceeds 37 degrees C or 100 degrees F.
- Do not use the device in bathroom or any humid places.
- If the device work disorder, please immediately cut power off, and contact the product after-service center.
- Do not use the device when you feel unwell.
- Confirm the device stop working when exchange the message heads.
- If the power wire is damaged, please ask the after-sales maintenance technician to change it.
- If the power cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

MAINTENANCE & STORAGE

1. To clean, wipe the machine with a soft cloth and neutral detergent only.
2. It is not allowed for any water or liquid into the device
3. Store it in dry place after use, avoid high temperatures, direct sunlight and corrosive gases.
4. Always unplug the unit if not in use







OPERATION

1. When inserted the plug in the socket, the power light is turned on.
2. Rotate the power switch / speed dial to right to start vibration massage and regulate vibrating intensity. Turn it to left to reduce intensity to power off.
3. Unplug the product after use or when they're not in use. Usage of

MASSAGE ATTACHMENTS


1. Flat massage head: Decompression massage, weight-loss massage for face, abdominal
2. Wavy massage head: Acupuncture points and deep tissue massage
3. Ball massage head: can be used with essential oil as you like
4. Micro filer massage head: Remove calluses and dry skin for feet. (Notice: Do not use mesh cover when micro filer massage head is attached.)
5. Mesh cover: Prevent hair involved in machines, protect hair safety.

DIRECTIONS FOR USE

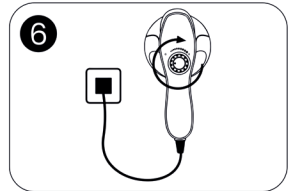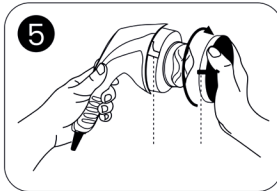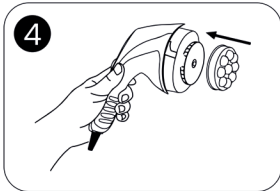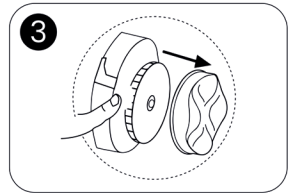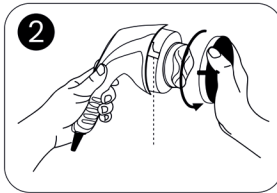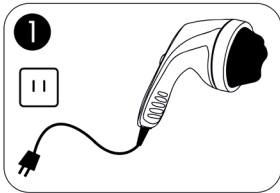1. Before you replace the massage head, remember to disconnect the power
2. Hold the host in your left hand, and pinch your right hand on the edge of the net. Twist to the inside, and the lock is loose and pulled out
3. Find the two small gaps between the head and the fuselage. Grab two gaps and pull out
4. Change the massage head you want and aim at the press
5. Hold the host in your left hand, and pinch your right hand on the edge of the net. The inside of the hood is pointed at the grip of the lock, and wiggle the card tightly
6. Plug in the power supply and turn on the switch. Just can use

PARTS AND FUNCTIONS

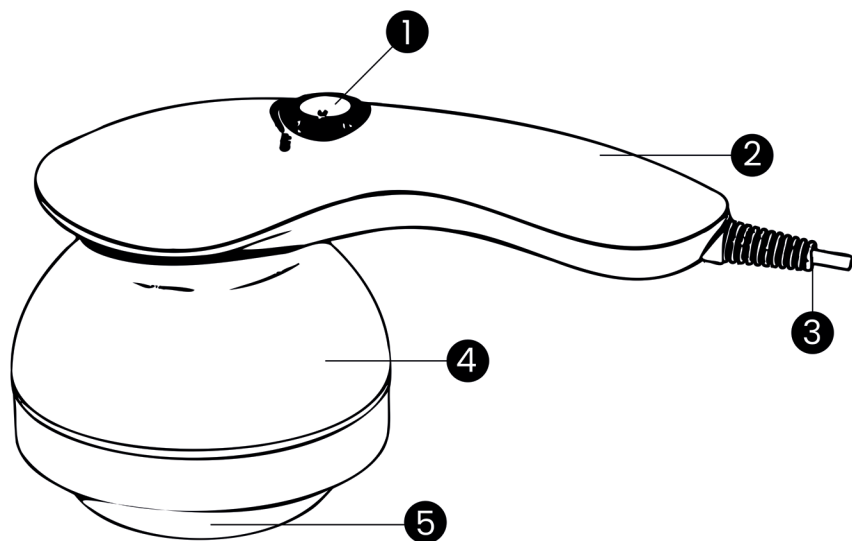1. Power Switch/Speed Dial
2. Hand Grip
3. Power cord
4. Main body
5. Massage Head

MULTIPLE MESSAGE ATTACHMENTS FOR CHOICE

TAIJI
MESSAGE
HEAD

BALL
MESSAGE
HEAD

ROLLER
MESSAGE
HEAD

FLEXIBLE
ANGLE
WHEEL


POINT ROLLER
MESSAGE HEAD

SMOOTH
MESSAGE
HEAD

SCRUB
MESSAGE
HEAD

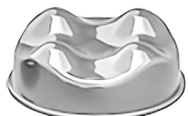WAVE
MESSAGE
HEAD

Our Body Massager is a new innovation that uses deep and powerful high frequency vibration and adjustable speed to help you smoothen, sculpt, massage and tone your body all from the comfort of your home.

Please read this section for the correct and safe use of this product before use. Here is a reminder of how to use products correctly to avoid danger or damage to anyone. After reading, please put the manual in a safe place so that you can check it at any time.

TECHNICAL DATA

SYMBOL	113860
INPUT	220-240V~ 50/60HZ
MAIN MACHINE MATERIAL	ABS
NUMBER OF ATTACHMENTS	8
SPEEDS LEVELS	SMOOTH ADJUSTMENT WITH A KNOB

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

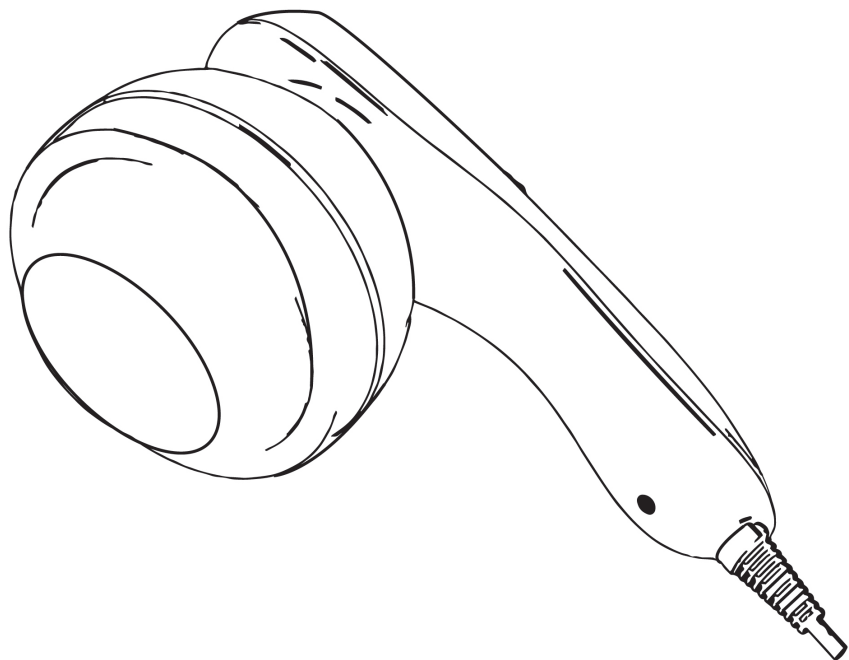
Malis B. Machoński Sp. K. declares that described in this manual **Maltec BODY MASSAGER SKINPULSE** meets all the requirements of the:

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