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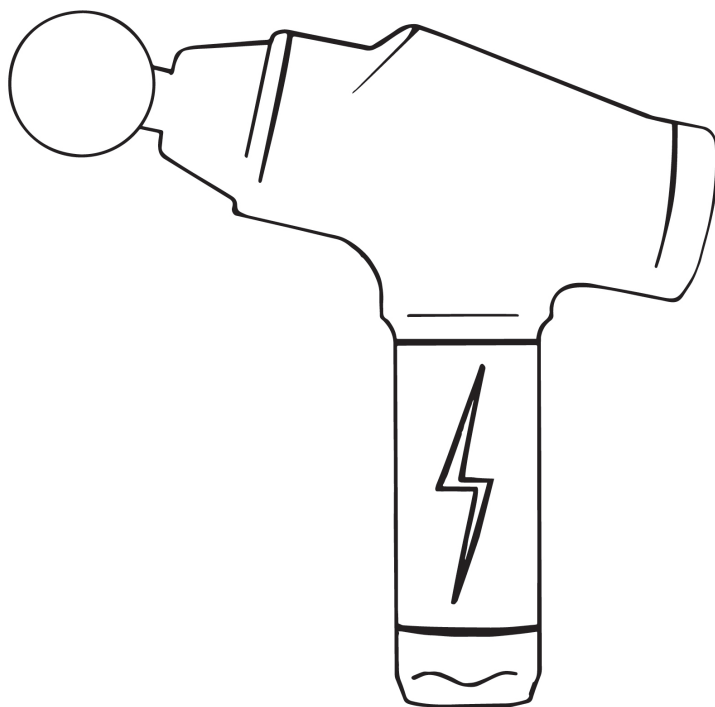
PL

PISTOLET DO MASAŻU
BODYBEAT FLEX

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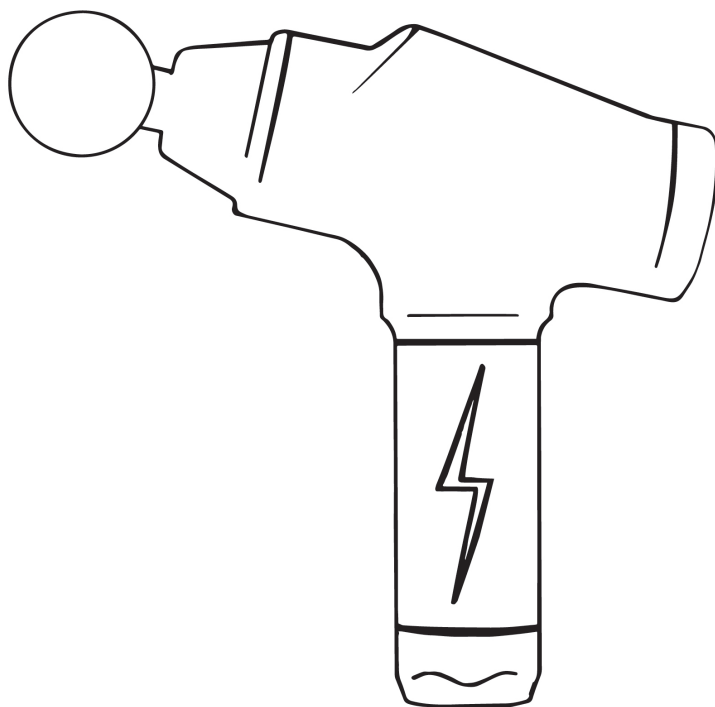
EN

MASSAGE GUN
BODYBEAT FLEX

**Before using this device, please read this user manual! Follow all safety instructions to avoid damage caused by improper use!
Keep this user manual for future reference.**

SAFETY INSTRUCTIONS

The massage gun is a device used for deep relaxation of muscles and fascia, suitable for physical massage after exercise and fitness. Thank you for purchasing our product. This manual will provide instructions for the effective and safe use of the product. Please read the instructions carefully before use and store them in a safe place for future reference. If you have ever passed on the product to others, this manual should be included.

USE TO

Suitable for physical massage after exercise and fitness

SAFETY PRECAUTIONS

Before use, it is necessary to understand the following content, otherwise it may cause damage or danger!

ABOUT THE USER GROUP AND LOCATION

- For use only by healthy adults other than the elderly;
- Only use the product on soft tissues of the body without causing pain or discomfort;
- Do not use on any hard or bony parts of the head or body;
- Whether under moderate pressure or under applied pressure, bruising may occur. Regularly check the massage area and stop using immediately if there are signs of pain or discomfort;
- Can only be used on dry, clean skin surfaces of body parts, gently pressing and moving on the skin;
- Keep fingers, hair, or other parts of the body away from the axis and back of the massage head, as they may be pinched;
- Do not use in areas with less muscle coverage near the bones, such as the tibia, dorsum of the foot, or dorsum of the hand;
- Do not use in or around sensitive areas such as eyes, teeth, breast augmentation, reproductive organs, prostheses, etc;
- Do not use the product in any possible fracture area;
- If the pain is persistent, immediately stop using and consult a doctor. This product is not a medical device and is not recom-

mended for medical use

- Not intended for medical use in hospitals;
- Before use, please check the host and connecting wires. If there is any damage, it must be replaced by the manufacturer or its designated professional personnel to avoid danger;

ABOUT INSTRUMENT OPERATION

- Use a massage head that can produce the most ideal effect. Do not place any objects in the ventilation holes of the product, do not immerse them in water, and do not let water enter the product
- Please prevent falling or improper use, which may cause permanent damage to the product. Before each use, the operating status of the host, battery, charger, etc. should be carefully checked to ensure that they can meet normal use.
- Do not change or modify the product in any way;
- Please do not let the product run or charge unattended, and do not start the device during the charging process
- Prohibit throwing equipment into fire;

- | | |
|--------------------------|-----------------------------|
| 1. Battery level display | 8. Replaceable massage head |
| 2. Decreased speed | 9. Handle |
| 3. Start/Stop | 10. control screen |
| 4. Pressure sensing lamp | 11. Battery level display |
| 5. 10 minute stop prompt | 12. Charging cable |
| 6. Gear display | 13. Charging port |
| 7. Speed increase | 14. Power switch |

- | | |
|---------------------|-----------------------|
| 1. Ball Head | 8. Arc-shaped head |
| 2. Flat Head | 9. Big feet head |
| 3. Fork Head | 10. Cold massage head |
| 4. Cylindrical head | 11. Small ingot head |
| 5. Air cushion head | 12. Thumb head |
| 6. Four claw | 13. Orange pot head |
| 7. Half round head | |

USAGE

Turn on the power switch and the machine will automatically enter 0 gear after starting up. Press the screen start button to enter 1 gear and start running. Press the screen+button to increase the gear, which will accelerate the speed and increase the force. Press the screen - button to decrease the gear, which will slow down the speed and decrease the force

SUSPEND

In working mode, press the Start/Stop key to change the machine to 0 gear and stop working. Press the Start/Stop key again to enter 1 gear and start working again

10 MINUTE SELF STOP

The working time for each operation is 10 minutes, and pauses during the operation are not counted as 10 minutes of working time. If the cumulative working time reaches 10 minutes, the equipment will automatically change to 0 gear and stop working. If you need to use it, please restart.

MEANING OF INDICATOR LIGHTS

After starting up, the default gear position is 0 and it is always on.

When charging, the gear indicator light changes to CH flashing. When fully charged, CH stops flashing. When the display shows LO, it indicates low battery. Please charge as soon as possible. After the battery is fully charged, the gear display will show CH and stop flashing. When charging, the machine cannot be started and used. When not in use, it should be charged at least once every 1-2 months. After being fully charged, it should be placed to avoid damage caused by battery depletion

PRESSURE SENSING LIGHT

1. When there is no pressing action, the indicator light goes out;
2. When the pressure is normal, the pressure sensing light turns green
3. When the pressure is too high, the pressure sensing light turns red
4. Please use this product when the pressure sensing light is green

BATTERY INDICATOR LIGHT

After turning on, the battery indicator light at the handle remains on, and the screen battery indicator light is 0-20-40-60-80-100. As the battery consumption increases, it gradually goes out. When the screen battery indicator light is only red, the battery level is below 20%. Please charge it in a timely manner

The following individuals are prohibited from using: implanting cardiac pacemakers, suffering from heart disease, strain back disease, pregnant women or shortly after delivery, patients with fractures, osteoporosis, joint disease, malignant tumors, skin trauma, or eczema. The use of equipment to replace or supplement medical treatment is also prohibited. If you have pain, are receiving medical treatment, or are using other medical devices, be sure to consult your physician before using the device. Individuals with physical discomfort should consult a doctor for advice before use

TECHNICAL DATA

SYMBOL	113857
POWER	45W
SPEEDS LEVELS	30
BATTERY	Li-Ion 2500mAh
CHARGING TIME	2-3h
USAGE TIME	3-4h
AMPLITUDE	12 mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

Li-ion

DECLARATION OF CONFORMITY

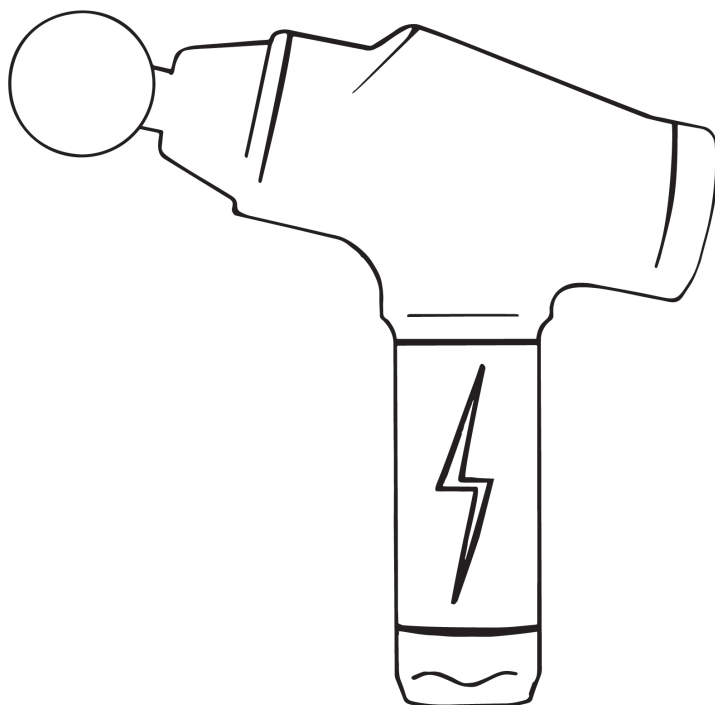
Malis B. Machoński Sp. K. declares that described in this manual Maltec MASSAGE GUN BODYBEAT FLEX meets all the requirements of the:

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