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
INSTRUCTION MANUAL EN

ULTRASONIC TOOTH SCALER
TOOTHIE X100

BEFORE USING THIS DEVICE, PLEASE READ THIS USER MANUAL! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THIS USER MANUAL FOR FUTURE REFERENCE.

PARTS OF THE PRODUCT

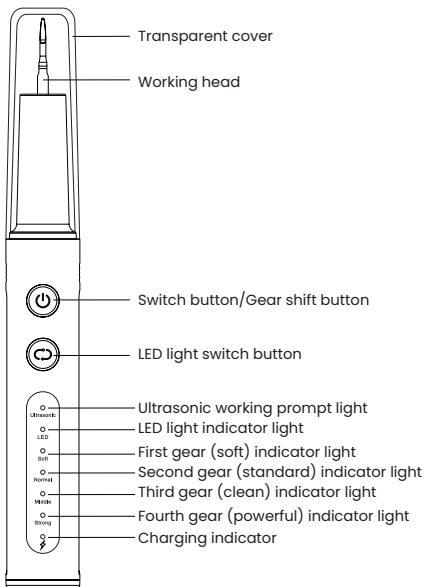SAFETY INSTRUCTIONS

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should be supervised to ensure they do not play with the device.
- Close supervision is necessary when any appliance is used by or near children.
- Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.
- Do not leave packing materials in reach of children because they can be dangerous.
- The product is not equipped with a power adapter. The customer's own adapter shall select output voltage of 5V DC.
- Before connecting the product to the power supply, check whether the voltage marked on the product is consistent with the local power supply voltage.
- If any component of the device is damaged, especially the power cord or the housing, stop using the device as this may pose a direct threat to health and life. It must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Do not immerse the charging port in water or any other liquid.
- The product cannot be used when charging.
- Please do not open the product or change the parts of the product without permission, so as

not to damage the product or cause electric shock, and the warranty after private demolition is invalid

- This product has no fire prevention function, please stay away from the fire source or the corrosive liquid.
- Patients with hemophilia and cardiac pacemakers are prohibited.
- No use by pregnant women and minors.
- No use by people with dental and porcelain teeth.
- No use by people with dental caries, periodontitis, gingivitis and tooth sensitivity.
- If the product falls into water during charging, please first unplug the power adapter from the power socket and then take out the product to avoid electric shock hazard.
- During charging, don't unplug the power plug and the charging plug with wet hands.



USER INSTRUCTIONS

Direction for use

When starting the tooth cleaning, do not aim the metal cleaning head at the sensitive parts of the human body, such as eyes, power supply, lamps and other dangerous equipment. Do not contact the skin with metal to avoid injury when opening the product.

1. Point the cleaning head at the top port of the sensor, align the bayonet with a wrench and tighten clockwise.
2. Replace or remove the cleaning head, align the wrench with the bayonet, and then release it counterclockwise to replace or remove the


cleaning head.

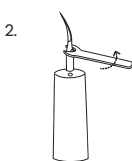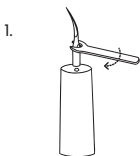
- Power on: Long press the power mode switch button to start, and the work indicator light is on.
- Switch mode: after turning on, start the last shutdown mode by default, touch the button again to switch the working mode, soft/normal/ middle/ premium mode in turn.
- LED light button: start the default LED light on, touch the LED light key to close, touch the LED light key again to open the LED light;
- Shutdown: press the power/mode switch button for 3 seconds, or automatically after 10 minutes.

Methods to clean your teeth


When cleaning the calculus, stains and plaque on the gums, please use a pointed tip to clean the calculus slowly and repeatedly.

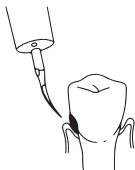
When removing large calculus and plaque from the neck and side of the tooth, please use a flat head for cleaning, and clean it slowly and repeatedly.

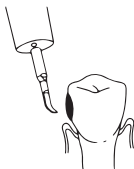
Longpress the one-button switch for above 2 seconds to turn on the power, and then according to the thickness of the tooth stone, press the power/mode switch button to adjust to the desired mode.

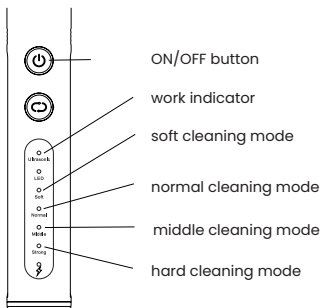
When cleaning dental stains and tartar, please use both sides of the cleaning head to avoid scratching the surface of the tooth; Place the cleaning head at an Angle parallel to the tooth surface or less than 15° , gently touch the lower part of the tooth stone, and move back and forth to shatter and shake off the calculus.

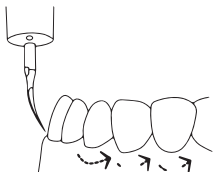Please use a very light force, will clean the head will move back and forth. Do not hold the cleaning head at one point to vibrate, otherwise it will cause damage to the surface of the tooth.

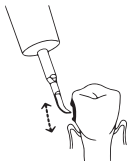Start from the place with little calculus, and use in the order shown in the figure above. Use an oral mirror for the unseen location.

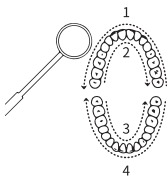CHARGING

NOTE: The product is not equipped with a power adapter. The customers's own adapter shall select output voltage of 5V DC.

NOTE: The product cannot be used when charging.

NOTE: If the charging port of the main units gets wet with washing water or water, please wipe it with a dry cloth before connecting to the charging cable and start charging after drying.

Connect the USB port at the bottom of the handle to the power adapter via the USC charging cable.

When charging, the charging indicator light flashes red. When fully charged, the charging light always lights on and stop flashing.

CARE AND MAINTENANCE

- Please wipe the shell with a soft cloth, do not use a corrosive detergent to wipe;
- Do not place this product in the sunlight, dust and polluted environment;

STORAGE

Always store the device:

- clean (see chapter CARE AND MAINTENANCE)
- dry
- protected against dust, frost, source of heat and damp
- out of the reach of children

Disconnect the plug from the outlet socket.

Give the machine time to cool off before you store the machine in an enclosed space.

Do not place other objects on top of the appliance.

TROUBLESHOOTING

fault phenomenon	The cause of the problem	Rx
The cleaning-head is not properly connected to the power shaft	There is a foreign body blockage in the cleaning head jack	Remove residual foreign body
	Power shaft adhesion to the foreign body	Remove foreign body
	The cleaninghead is not genuine or mismatched	Buy genuine cleaningheads
	The direction is not correct	Position it correctly as indicated
Unable to start	Too low or run out	charging interval
	be responsible for	Pull off the charging line

Cleaning head vibration is weak	Cleaning head is connected in correctly or improperly connected	Correct installation
	There are foreign matter between the cleaning head and the power shaft	Remove foreign objects
Clean the head without shaking	Untriggered on key	Press the power button again
	The operation head is not tightened	Tighten the operator head with a wrench
	Run for over 10 minutes	Press the power button again
	Battery power is insufficient	Please charge it up before you use it
The auxiliary lamp does not work	The LED button is not triggered	The LED button is not triggered. Press the LED button again
Charging status is not displayed after connecting to the charger	This plug is not well connected to the socket	Check that the power supply is switched on
No charging display	Power out or low	Normal charging can be restored

TECHNICAL DATA

INPUT POWER	3.5W
INPUT VOLTAGE	5V DC, 1A
BATTERY	1200mAh, 3.7V
OPERATING FREQUENCY	36000Hz
HEAD LIGHT	yes
DISPLAY	yes
METAL HOUSING	yes
4 OPERATING MODES	yes
IPX7 WATERPROOF	yes
INCLUDES	4x heads, 4x auxiliary tools, charging cable, head cover

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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