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
INSTRUCTION MANUAL EN

WATER FLOSSER
TOOTHIE X70

BEFORE USING THIS DEVICE, PLEASE READ THIS USER MANUAL! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THIS USER MANUAL FOR FUTURE REFERENCE.

PARTS OF THE PRODUCT

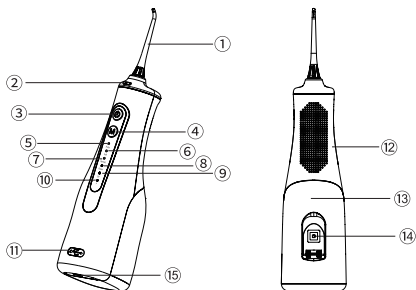- 1.Nozzle
- 2. Nozzle release key
- 3. Power button
- 4. Mode button
- 5. Pulse mode indicator
- 6.Soft mode indicator
- 7. Normal mode indicator
- 8. Power mode indicator
- 9. DIY mode indicator
- 10. Power indicator

- 11. Type-C charging port
- 12. Main body
- 13. Water tank
- 14. Disinfection lamp
- 15. Charging port

SAFETY INSTRUCTIONS

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should be supervised to ensure they do not play with the device.
- Close supervision is necessary when any appliance is used by or near children.
- Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.
- Do not leave packing materials in reach of children because they can be dangerous.
- The product is not equipped with a power adapter. The customers's own adapter shall select output voltage of 5V DC.
- Before connecting the product to the power supply, check whether the voltage marked on the product is consistent with the local power supply voltage.
- If any component of the device is damaged, especially the power cord or the housing, stop using the device as this may pose a direct threat to health and life. It must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Do not immerse the charging port in water or any other liquid.
- The product cannot be used when charging.
- Please do not open the product or change the parts of the product without permission, so as


not to damage the product or cause electric shock, and the warranty after private demolition is invalid.

- This product has no fire prevention function, please stay away from the fire source or the corrosive liquid.
- Don't fill the container with water above 40°C to clean your mouth to avoid scorching.
- If the product falls into water during charging, please first unplug the power adapter from the power socket and then take out the product to avoid electric shock hazard.
- Don't throw the product into the fire or heat it.
- Don't charge the product in a hot or humid place.
- Don't place the product in these places for a long time and meantime, avoid direct sunlight.
- During charging, don't unplug the power plug and the charging plug with wet hands.
- If the product is not used for a long time, please unplug the power plug and the charging plug to avoid fire or other hazards.
- Do not stare directly at the UV lamp.

USER INSTRUCTIONS

Nozzle mounting


Insert the nozzle into the nozzle holder on the top of the main unit until it clicks.

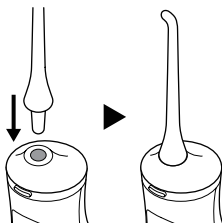Injection of cleaning solution

Open the lid of the water tank filling port, inject clean water into the tank and then cover the water filling port lid.

If you are using a concentrated type of mouth wash, then please dilute it with water and fill the mouth wash.

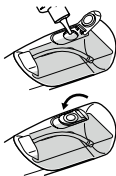Usable cleaning solution


On the right is a list of liquids that can be put into the can and used for cleaning. Do not use any other liquids as this may cause malfunction.

- tapwater
- warm water (40°C or lower)
- Mouthwash

Mode settings

Press the mode switch button on the front of the main unit and select the desired water pressure. Every time it is pressed. It switches in the order of Pulse-→Soft-→Normal-→Power-→DIY.

After turning on the device using the power button, press and hold the mode button for 3 seconds in DIY mode. The indicator light will flash, entering the DIY settings mode. Press the mode button to cycle through 10 settings. Once the desired setting is selected, the light will flash 5 times to confirm completion.

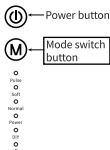Power on

Insert the lip of the nozzle into the mouth and press the power button. The indicator lights up to indicate that cleaning has started.

Cleaning


By leaning forward and opening your mouth slightly during cleaning, you can flush cleaning water from your mouth.

*Please note that if you open your mouth wide and use ii, the washing water will splatter around.

The main body is equipped with a waterproof mechanism, so you can use it with confidence even if it gets wet with washing water.

The direction of the nozzle can be freely rotated 360°. Please use ii by adjusting it to the direction that is easy to use.

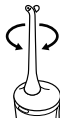If the lip of the water supply hose inside the tank rises above the water surface, washing water will not come out. Adjust the direction of the main body and use it with the tip of the hose submerged.

Point of care

Aim for the water between the teeth and between the gums. then move slowly along the gums from the molars (back te-

eth)to the front teeth.Clean the surfaces and backs of teeth.


Stop

When you want to stop the operation of the main unit. such as when the cleaning water is used up, press the power button again.

How to remove the nozzle

Pull the nozzle up while pressing the nozzle release button.

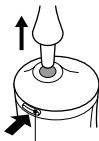CHARGING

NOTE: The product is not equipped with a power adapter. The customers's own adapter shall select output voltage of 5V DC.

NOTE: The product cannot be used when charging.

1. Connect the USB-C port of the supplied charging cable please remove the cap of the charging port of the main unit and connect it to the charging port. Connect the other end (USB-A) connector of the charging cable to your USB output charging adapter.

If the charging port of the main unit gets wet with washing water or water, please wipe it with a dry cloth before connecting the charging cable and start charging after drying.

2. When charging begins, the power indicator flashes continuously to show that charging is in

progress.


3. When charging is complete, the power indicator stays lit.

CARE AND MAINTENANCE

- After use, open the tank and discard the remaining water in the tank each time. After that, press the power button to activate this product, and drain the remaining water from the nozzle. When water stops coming out, press the power button to stop using.

- If dirt adheres to the inside of the tank, put a small amount of lukewarm water of 40°C or less in the tank and wash it a little, then perform the above.


- To keep it clean, rinse the nozzle with water after each use and let it dry.

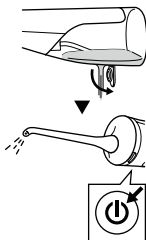If no water sprays out for a long time. Please I unplug the nozzle and press the water valve unplug the nozzle and press the water valve in the water outlet intermittently With a tool With a diameter of less than 6mm until water sprays out.

There is an air inlet in this part of the water tank, and when filled with water and turned upside down, there may be a small amount of water flow that is normal. When air enters during startup, ensure that the air pressure inside the water tank is balanced.

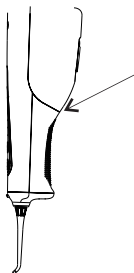STORAGE

Always store the device:

- clean (see chapter CARE AND MAINTENANCE)
- dry
- protected against dust, frost, source of heat and damp
- out of the reach of children

Disconnect the plug from the outlet socket.

Give the machine time to cool off before you store the machine in an enclosed space.

Do not place other objects on top of the appliance.

TECHNICAL DATA

POWER	6W
VOLTAGE	5V
PRESSURE	2.8-9.0bar
TANK	300ml
BATTERY	2000mAh, 3.7V
5 OPERATING MODES	yes
UV STERILIZATION	yes
IPX7 WATERPROOF	yes
INCLUDES	6 interchangeable nozzles, charging cable, case

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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