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MASAŻER DO RĄK
HANDFORCE WHITE

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INSTRUCTION MANUAL EN

HAND MASSAGER
HANDFORCE WHITE

PLEASE READ THIS USER MANUAL BEFORE USING THE DEVICE! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THE USER MANUAL FOR FUTURE REFERENCE.

SAFETY INSTRUCTIONS


Please follow the doctor's advice for the following patients:

The patients who have had hand surgery, swelling and inflammation, congestion, severe arthritis, trauma and other hand diseases;
People with heart, brain and lung diseases;
Blood diseases with bleeding tendency;
Osteoporosis and soft tissue diseases, stroke;
Infants and pregnant women.

1. Use Type-C plug connection.(Note the shape of the insertion port)
2. Before use, Please make sure that you have no jewelry(watch ring, Bracelet ring, etc.) in your hands, so as to avoid sharp objects piercing the airbag or squeezing the airbag and causing jewelry damage.
3. It is forbidden to use for a long time. Normally, massage time should not exceed 30 minutes.
4. Please do not use this product while driving.
5. Please do not use it after getting drunk or taking sleeping pills.
6. Please use DC5V1~2.5A charger.
7. This product is specially designed as household electronic health care products. It is not suitable for diagnosis, treatment and other professional or commercial purposes.
8. Please keep this product away from water, fire and corrosive environment.
9. Please do not put sharp objects into this product.
- 10 Please do not use when hands are wet.
11. If the power cord is damaged, stop using it immediately. Contact the product service center for maintenance. Do not disassemble the product at will.
12. If you feel unwell, please stop using it in time.
13. Please do not wrap products with towels, quilts, clothes, etc. when you use it.

PRODUCT DESCRIPTION

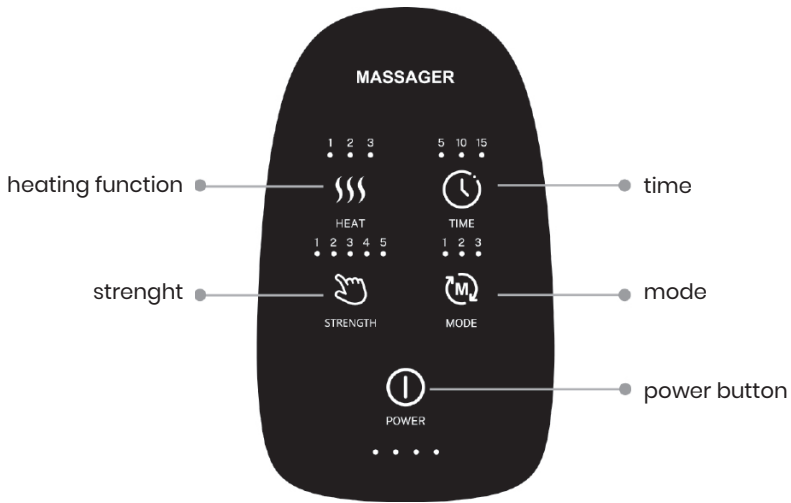USER INSTRUCTIONS

1. Stretch the back of your hand upward into the machine.
2. Place the palm of your hand on the circular raised area in the middle of the massage pad.
3. Press the power start button on the panel, and all button indicators will light up.
4. Press the strength button on the panel to select different strength modes according to different sizes of hands.
5. Press the mode button M to select gentle mode, soothing mode, and strong mode from low to high, a total of 3 different strengths and different massage modes.
6. If you need to start the heating function, press SS5 to turn on the low heating range, then press the medium range, then press the high range, and press again to turn off the heating.
7. Press the time button to select 5 minutes, 10 minutes, 15 minutes, different time modes.
8. If you need to stop using or turn off the power, please press I.

- Try not to use it while charging.
- When the product is charging, the fingerprint logo flashes until it is fully charged and stays on.
- Apply hand cream fully to your hands before using this device, which can help improve hand lines and skin.
- Do not use it in heating mode for a long time to prevent low-temperature burns.

MAINTENANCE AND CLEANING / STORAGE

1. Use towel to remove dust and other dirt, and place it in a dry and ventilated place.
2. Do not place in high temperature, humidity Maintenance Soft tissue diseases; stroke; infants and pregnant women; and storage and direct sunlight.
3. When the storage period exceeds 6 months, it should be checked out and power-on.
4. When plugged into the USB socket, the machine can be used. But we strongly recommend against it.

TROUBLESHOOTING

If there is a malfunction during the use of this product, please check it as follows. If the fault can not be eliminated or the product still cannot run, please turn off the power and contact the product service center for maintenance.

Basic Faults	Possible causes	Troubleshooting Measures
The product can't start working	Check the remaining power of the battery The power may be weak.	Charging
The product cannot be re-charged	Is the power adapter connected to the data line correctly?	Re-plug until LCD power display flickers
	Is the data line and USB socket intact?	Check the socket If there is any damage, please contact the manufacturer
	Are the specifications of adapters and data lines in line with operational requirements?	View specifications and proofread the technical parameters provided in the instructions
The product can be turned on, but can't work	The remaining power may be weak.	Charging

TECHNICAL DATA

RATED VOLTAGE	DC 5V
BATTERY CAPACITY	2000mAh
RATED POWER	7.5 W
PROTECTION TIME	15 min
CHARGING	5V 1A
WEIGHT	0.75kg
MATERIAL	Antibacterial and deodorizing fabric

POWER SUPPLY NAME	Lithium battery
OUTER SIZE	250x154x111 mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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