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PL

MASAŻER DO STÓP
FOOTPULSE BLACK

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INSTRUCTION MANUAL EN

FOOT MASSAGER
FOOTPULSE BLACK

PLEASE READ THIS USER MANUAL BEFORE USING THE DEVICE! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THE USER MANUAL FOR FUTURE REFERENCE.

SAFETY INSTRUCTIONS

In order to avoid Improper use, resulting in burns, electric shock and other dangerous situations. When using this product, please be sure to pay attention to the following matters.

1. After use and before cleaning the product, please remove the power plug from the socket.
2. If the product falls into the water, do not put your hand into the water to retrieve the product, please unplug it from the socket first.
3. Do not use this product in the bath or shower.
4. Do not place or store the product in any place that will cause the product to be immersed in water.
5. Do not place or immerse the product in water or other liquids
6. Do not use sharp objects on this product.
7. Please keep the product dry and do not act in a humid environment.
8. Please unplug the product from the socket when no one is watching.
9. Please do not let children play with this product.
10. Do not use the product if the power cord or plug is damaged, not working properly or has fallen into the water.
11. The power cord must be kept at a safe distance from the heating surface. Please do not place the product near the source of fire and flammable places to avoid failure or danger.
12. This product has a refreshing effect, please do not use this product when sleeping.
13. Do not insert or drop objects into any openings of the product.
14. Do not use spray products or places filled with flammable bodies to operate this product.
15. Do not place the product under a blanket or pillow to avoid overheating the product and causing a fire, short circuiting that wires or injuring people.
16. Do not drag the power cord to move the product, or use the power cord as a grip.
17. When unplugging the power, please make sure the product has been turned off first to avoid damage to the product.
18. This product is not for any medical use.
19. Do not subject this product to strong impact.
20. Do not use parts not provided by the manufacturer. Be sure to use the product in accordance with the instructions described in the user manual.

If the power cord is damaged it must be replaced by the manufacturer or a professional repairer.


If you have the following conditions, please consult a doctor before using this product:

- During pregnancy, during menstruation or shortly after parting of the night
- Those who use medical electronic devices such as cardiac pacemakers that are implanted in the body
- People with blood circulation disorders or thrombosis
- Have osteoporosis
- Skin infections
- High fever

- Traumatic injuries or superficial diseases
- Doctor's orders to recuperate or feel unwell
- The surface of the product will be hot during use. so people who are not sensitive to heat must be careful when using the product.

1. Each use time of no more than 15 minutes (the product automatically turns off after completing 15 minutes of massage program), the use of too long may lead to overheating, thereby shortening the product life, such as the product overheating, please immediately stop using, first cool the product before the next operation.
2. In order to reduce the impact of the hedge current on the product, please do not frequently switch the power supply within a short period of time to extend the life of the product.
3. Please use the power adapter configured for this product, using other power adapters may cause damage.
4. Do not use the damaged power supply or tamper with the circuit of this product.
5. If a power failure should be unplugged to avoid the sudden restoration of power to cause injury.
6. If you feel physical abnormality during the use of this product, please stop using it immediately and consult a doctor.
7. This product is not recommended for people with physical disabilities, sensory or neurological deficiencies or lack of experience and common sense (including children, the mentally ill, the physically paralyzed, etc.) to use, if necessary, must be used under the supervision and guidance of their guardians.

PRODUCT DESCRIPTION

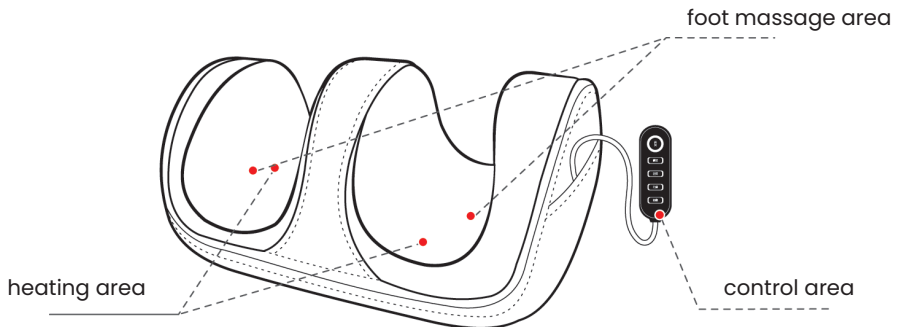
USER INSTRUCTIONS

Prepare before use

First check whether the adapter power plug is intact, and then well placed the product. Making massage aligned with the using part, then make DC plug on the right side of the product is docked with the adapter, and then plug adapter into the socket.

Remote control description

Button	Button Description
Power	Short press to turn on. Short press again to turn off. The red light is turned on by default, and the medium speed is turned on.
Mode	Default automatic mode, short press mode one, then press mode two, short press mode three for the third time. And so on.
Reversion	Default positive rotation. Short press reverses the rotation, the second short press restores the positive rotation.
Intensity	Default medium, short press high, short press low again, short press medium for the third time.
Heating	The heating function is turned on by default, short press to turn off the heating, short press again to turn on the heating.
• It is recommended to massage the same part of the body for not more than 15 minutes each time, if you feel uncomfortable during the massage, please stop using it immediately. • When the massage pillow is not in use, make sure that the power plug is unplugged from the outlet. • In the event of a power failure, turn off the main switch and unplug the power supply to prevent the product from being re-operated when the power supply is suddenly restored.	

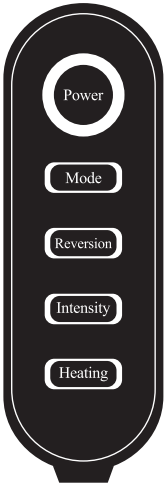MAINTENANCE AND CLEANING

Maintenance method

1. Usually disconnect the power when not in use and place the product in a safe, dry and cool place.
2. Avoid contact with sharp objects to avoid puncturing or damaging the appearance of the product.
3. Please do not store the product in high temperature, near flame or direct sunlight, otherwise it will cause fading and damage.
4. Avoid contacting the product and accessories with benzene, thinner or any other corrosive liquid.
5. When storing the product, do not tie the product with any power cord or wire or hang the product with wire.

Cleaning method

1. Before cleaning, please unplug the power supply and wipe gently with a soft, slightly damp sponge or soft cloth.
2. The product must not be immersed in any liquid, not directly rinse with water or any liquid directly.
3. Do not use sharp brushes, gasoline, kerosene, furniture cleaning fluid, etc. to wipe the product.

TROUBLESHOOTING

Failure phenomenon	Solution
Cannot start the massager	Please confirm that the power plug is corectly inserted into the socket. Please confirm that the socket switch has been turned on. Please confirm that the massager switch button has been pressed.
Massage suddenly stop in progress	Please check whether the power plug Is unplugged from the socket or loose at the power adapter The 15 minutes set by the outomatic timing device has arrived. Overheating protection device to start, please turn off the main power supply, and wait for the massager to cool down.

TECHNICAL DATA

POWER	24W
INPUT (AC ADAPTER)	AC 110-240V / 50-60 Hz
INPUT (CAR ADAPTER)	DC 12 V
RATED TIME	15 min
SIZE	365x195x255mm
WEIGHT	1.8kg

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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