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PL

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INSTRUCTION MANUAL EN

CUPPING MASSAGER
BODYSCULPT

PLEASE READ THIS USER MANUAL BEFORE USING THE DEVICE! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THE USER MANUAL FOR FUTURE REFERENCE.

SAFETY INSTRUCTIONS

1. Those with bleeding or bleeding tendency and varicose veins are prohibited.
2. Pregnancy, menstruation, skin damage, skin sensitivity, postoperative patients, patients with high fever or critical diseases are prohibited Use.
3. Bleeding prone diseases, such as thrombocytopenia, leukemia, hemophilia, Henoch Schon lein purpura and other diseasesIt is forbidden for patients with blood diseases.
4. Excessive fatigue, drunkenness and staying up too late are prohibited (2 hours before and after meals are also prohibited).
5. The 12 gear red light heating model has a heating surface, which is not sensitive to heat, and should not be used by infants or people who can't take care of themselves.
6. When using, it should be scraped in one direction, not one up and one down or one left and one right.
7. When in use, those with poor local microcirculation will have Sha, and the „Sha marks“ will be absorbed by themselves, which is a normal situation. For daily use, it is mainly used for gentle massage. It is not required to get rid of Sha. It can be scraped until the pores are open, and the local skin is flushed and hot.
8. The use time shall not be too long, and continuous large-area massage is not allowed.
9. After use, the blood vessels on the skin surface will dilate. You should rest for 45 minutes before bathing.

Note: This product is not a medical device and cannot represent that the medical device can achieve the therapeutic function.

PRODUCT DESCRIPTION


USER INSTRUCTIONS

After starting the product, select the product function according to the needs, and start using it after applying essential oil to the required parts of the body (For sensitive skin, please test on the inner side of the arm or the root behind the ear before using the essential oil. If there is redness, swelling and itching, stop immediately use).

When scraping, apply some essential oil to the required parts (do not add more essential oil, let alone add essential oil to the machine), Turn the negative pressure on a small gear and pull it in one direction after sucking the skin. The scraping and swabbing of each part is generally less than 20 times, with the strength from light to lightHeavy, from small to large.

Note: the machine is equipped with large filter cotton, which can effectively prevent essential oil from entering the machine and prolong and protect the service life of the machine. each After using the essential oil for the first time, please replace the filter cotton in time. The filter cotton can be recycled after cleaning and air drying. Negative pressure during cupping Turn to the appropriate gear, suck the skin for 3-5 minutes, tilt the machine and pull it up or press the pressure relief key to easily pull it out and cupping No essential oil is required.

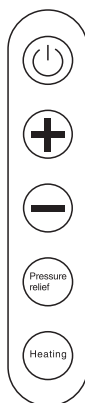This machine is a charging model: it is operated without a line after being fully charged, and cannot be started during charging The white light flashes when charging, and lights up when it's full.

This machine is equipped with a 20 minute timing function. The digital display screen displays a 20 minute countdown and automatically turns off the power after 20 minutes Source.

When the power switch is set up, press and hold for 2 seconds to start the machine, the digital 00 on the digital display screen will light up, and the 20minute countdown will be displayed after 2 econds; Press and hold for 2 seconds to turn it off, and the white light on the digital display screen goes out.

[+][-] gear strengthen/weaken keys, short press to adjust the negative pressure gear, there are gears 0-12 (gear 0 is no negative pressure), starting up machine defaults to 3 gears. Each time it is adjusted, the digital display screen displays the corresponding gears, and the countdown of 20 minutes will be displayed again after 2 seconds.

[pressure relief] it is effective when the negative pressure is turned on. Press it briefly and unload the gas for 1 second, which can easily solve the problem of difficult tank lifting.

[heating] short press to adjust the temperature. The heating function (i.e. gear 0) is not turned on by default when starting up. There are gears 0-12. Each gear is added, When the temperature is increased by one degree, the red light of the chassis lights up when the heating function is turned on, and the temperature is about 45 degrees, which is shorter when it reaches the 12th gear Press to return to gear 0, the digital display screen will display the temperature gear, and the 20 minute countdown will be displayed again after 2 seconds.

The charging time and service time of the product are determined according to the size and function of the user's charger.

In the above operation states, short press the power key to suspend the machine and enter the standby state, and short press again to restore the working state.

MAINTENANCE AND CLEANING

1. Please put the product in a safe place to avoid falling. Do not put heavy objects on the product to prevent damage.
2. Please store in a dry and well ventilated place to avoid strong sunlight, corrosive gas

and liquid erosion.

- 3. Please keep this product out of the reach of children.
- 4. When not in use for a long time, turn off the power supply and keep it fully charged every two months.
- 5. In case of technical failure, please contact technicians for maintenance in time. Do not disassemble the product without permission.
- 6. Before cleaning, the power supply must be cut off and the cleaning can be started only after the machine is cooled.
- 7. Clean with soft and slightly wet sponge or soft cloth. Do not clean the product with corrosive chemical cleaner.
- 8. **This machine is not waterproof. It is strictly prohibited to use in humid environment, and can not be washed with water to avoid water or other liquidsThe body enters the machine body to avoid short circuit or electric shock during use.**

TECHNICAL DATA

POWER	8W
CHARGING	5V 2A
VOLTAGE	DC 3.7V
MATERIAL	ABS PC
SIZE	200 x 90 x 90 mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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