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MalTec Health Massage Gun
PRO-GUN3800 7w1
Instruction Manual

SPECIFICATIONS

8. Powerful Brush-less High-Torque Motor with Quiet Glide Technology
9. Removable Battery Release
10. Rechargeable Lithium-Ion Battery (2+ Hrs of Use Per Charge)
11. Multi-Speed Settings for Optimal Pressure
12. Cordless and Lightweight "Includes AC wall charger
13. ON/OFF Power Switch
14. Interchangeable Applicators (x4)

USAGE INSTRUCTIONS

Please read prior to use

Charging

- Fully charge battery for up to six hours before first use either attached or detached from the device.
- To charge, connect the dC end of the supplied 12V adapter to the battery's charging port and plug ac side cable into wall outlet.
- LED band on battery lights pulsate to show battery level and indicate active charging.
- With battery attached, LED lights on unit flash to show battery level and indicate active charging.
- Full charge is indicated when all LED lights stop flashing and remain fully illuminated
- Battery may be recharged at any time and at any battery level.
- It is not recommended to fully deplete the battery to the red led level.
- The average run time is 2+ hours depending on speed level and pressure applied during use.

Operating

- Align and attach battery into base of unit-ensure it is securely locked in place.
- Power device on by moving the power switch on the bottom of the battery to the on position power device off by moving the power switch to the OFF position.
- With power switch in ON position, press the power button once to turn device on to level one twice to turn device to level two and three times for level three. Pressing the power button again will turn the device off. The blue LED lights will indicate the corresponding power level.
- With unit off, fully insert applicator by pressing in. Remove by pulling it out.
- Massage desired body part at the desired speed applying pressure as needed.
- To remove the battery, lightly press the battery release button with thumb. Firmly pull downwards on the battery handle with the other hand.


Maintenance, cleaning storage

- Use a slightly damp towel to clean device and a dry with a soft cloth.
- For storage or travel, remove battery and store in supplied carton box.

MASSAGE HEAD REPLACEMENT

if you want to remove the battery, gently press the battery release button with your thumb as shown in the picture, and pull the battery/handle with your other hand

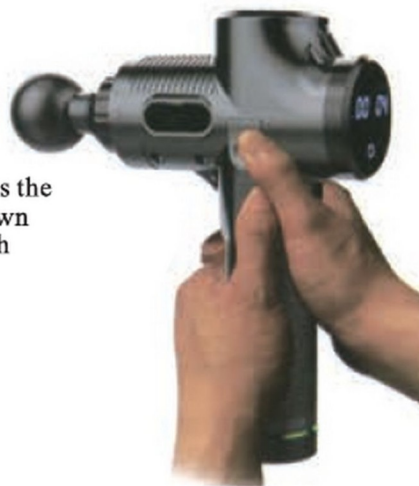DISCLAIMER

Warnings & Caution

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 30 minutes before re-using.
- Unplug the charger from device battery after charging or prior to use.
- For adult use only do not use if injured. Consult your doctor before using this product.

Safety Instructions

Danger: To reduce risks of electric shock, fire, and personal injury, this product must be used in accordance with the following instructions.

- For adult use only.
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the Massage Gun only on the body soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
- Use only the applicators that provide the best desired effect.
- Bruising can occur regardless of control setting or pressure applied Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the Massage Gun.
- Do not immerse in water or allow water to enter the ventilation holes of th Massage Gun.
- Do not drop or misuse the Massage Gun.
- Only recharge with supplied 26V DC Massage Gun charger.
- Carefully examine device and battery before each use.
- Do not tamper or alter Massage Gun in any way.
- Never leave the Massage Gun operating or charging unattended.

DISCLAIMER

Please do not use the massage gun or any vibration device without first obtaining approval from your doctor if any of the following apply.

Pregnancy diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery epilepsy or migraines, herniated disks, spondylolisthesis, spondylol-ysis, or spondylosis, recent joint replacements or IUDs, metal pins or plates or any concerns about your physical health frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in number of cases. This must be done on the advice and in accompaniment with a doctor specialist or physiotherapis.

CE declaration

Malis B. Machoński sp. K. Declares that the MalTec Health PRO-GUN3800 7in1 massage gun meets all the requirements of the EC Directive 2014/30 / EU and meets the following standards:

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