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PL

WODNY MASAŻER DO STÓP PRO-HYDRO500WS

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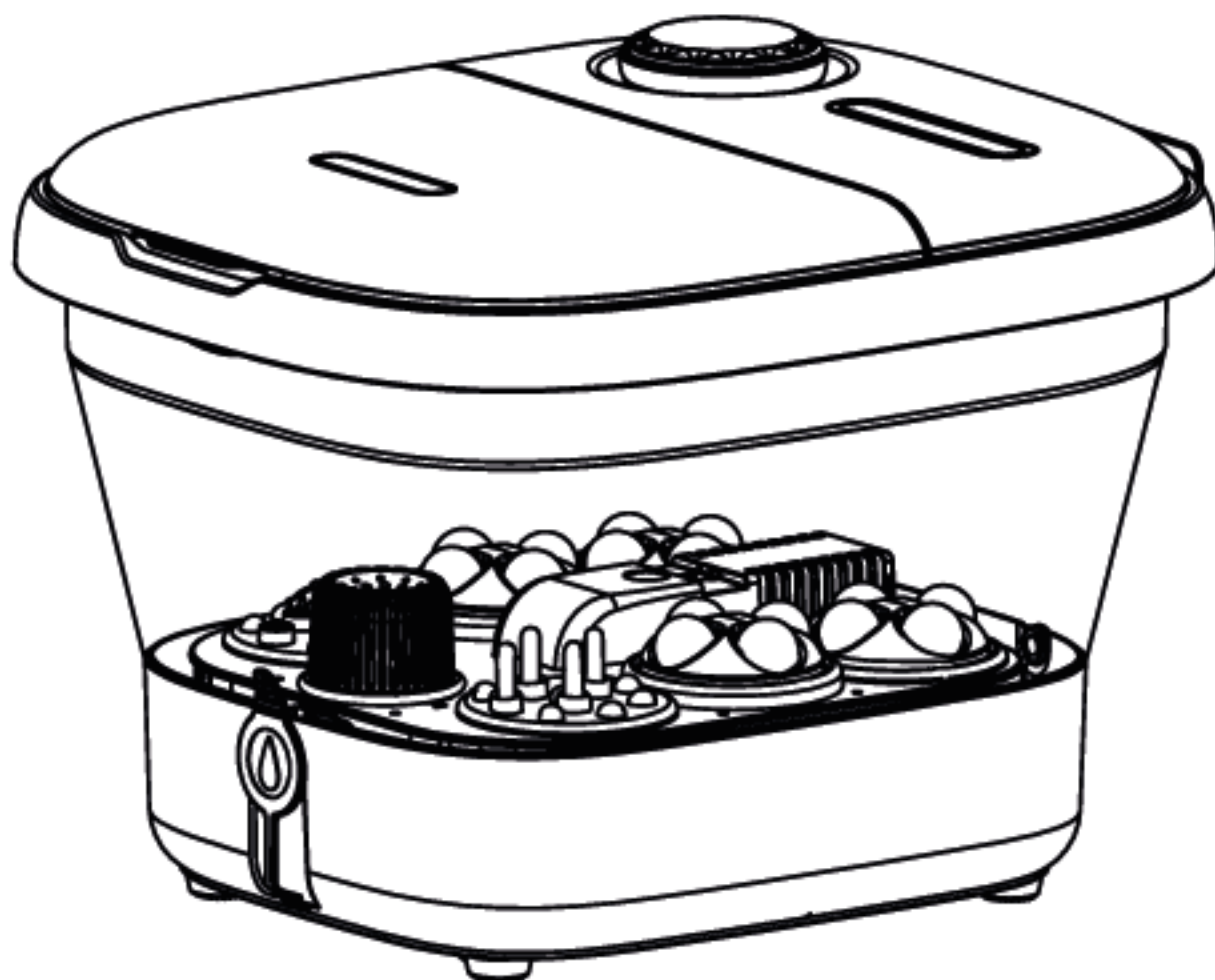
ENG

WATER FOOT MASSAGER PRO-HYDRO500WS


Thank you for your purchase. Please carefully read the manual carefully before, especially the precautions, so as to use this machine correctly. Please keep this manual properly so that you can read it at any time when necessary.

**Please use it according to the doctor's advice for the following people**


(Please listen to your doctor's advice and instructions before using any of the following conditions.)

- Psychiatric patients
- Implanted pacemaker or artificial organ in vivo
- People with cardiac dysfunction and neurogenic hypertension
- Patients with uncured cerebral haemorrhage and severe thrombosis
- Patients with haemorrhagic diseases, septicemia, fever, etc.
- Women during menstruation or pregnancy
- Those who are blunted or unconscious in response to temperature
- Those with inflammation, trauma or skin burns on the feet
- Children and the mentally handicapped must be used under the supervision of an adult

**PRODUCT PARTS**

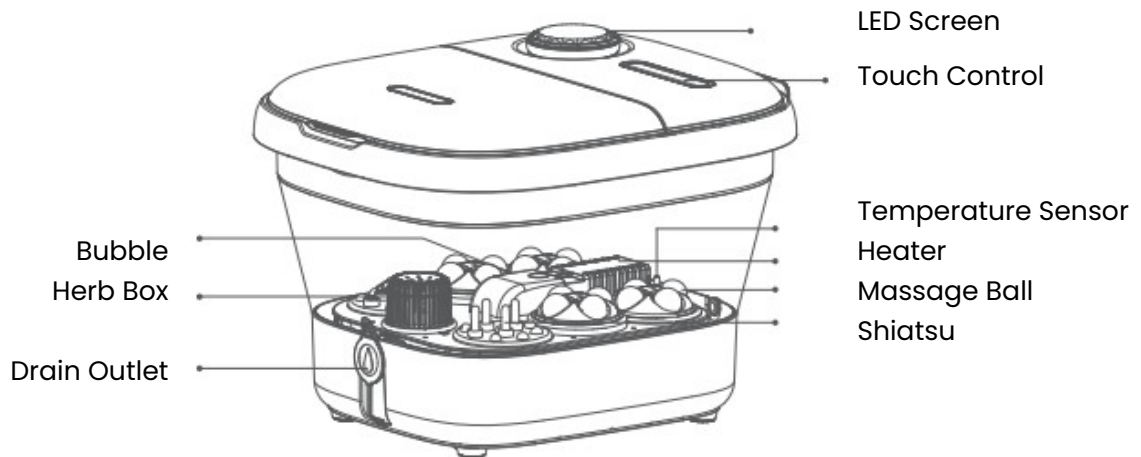Unfold

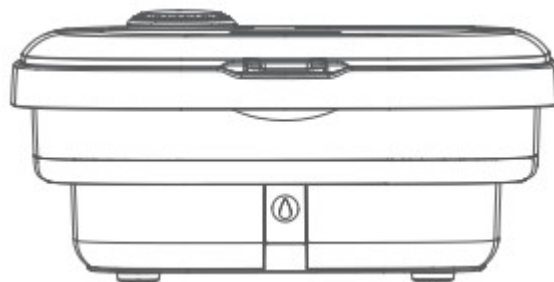Fold

## UNPACKING PARTS

1 Machine

1 Manual

## FUNCTION BUTTON

Control the whole machine on/off;

Control water temperature;

Control massage on/off (Shower head meanwhile on/off);

Control bubble on/off;

Control timing up to 60 minutes, 10 minutes increments;

**Attention:** Please unplug the product when not in use for a long time and store it in a dry place

## OPERATION INSTRUCTION

### PLUG IN THE POWER SUPPLY

- Touch "POWER", the prompt tone appears, start heating. The screen shows the actual water temperature. The default target temperature is 42°C and the working time is 60 minutes;
- Touch "TEMP", adjust the temperature from 35°C-48°C, and display the actual water temperature after setting for 3s;
- Touch "MASSAGE", turn on continuous massage. Touch "MASSAGE" again, turn on intermittent massage. Touch "MASSAGE" again, turn off the massage;

- Touch "BUBBLE", turn on the bubble. Touch the "BUBBLE" again, turn off the bubble;
- Touch "TIMING", display the remaining time of foot spa. Touch "TIMING" again to set your desired foot spa time, increasing each time by 10 minutes, up to 60 minutes;

## PRODUCT CARE AND STORAGE

- Before cleaning, please confirm that the power switch is turned off and the power plug is unplugged from the socket.
- Please wipe the surface of the product with a soft cloth.
- Clean the dirty area with a soft cloth dipped in water or a neutral solvent, and then wipe it dry with a dry cloth.
- Make sure the power switch is off and the power plug is removed from the socket
- Keep the product clean and away from places with high humidity or heat and volatile gases.
- If the products not used for along time, please put it in the packaging box to prevent dust from falling into it.
- Do not remove and insert the power plug with wet hands.
- Do not rinse the product with water or spill water or other liquids inside the product
- Do not use corrosive cleaning agents such as benzene or thinner to clean the product, otherwise it may cause the product color to fade and the product structure to be corroded and cracked.
- Do not place the product in direct sunlight or high temperature to avoid color fading.

## UNFOLD

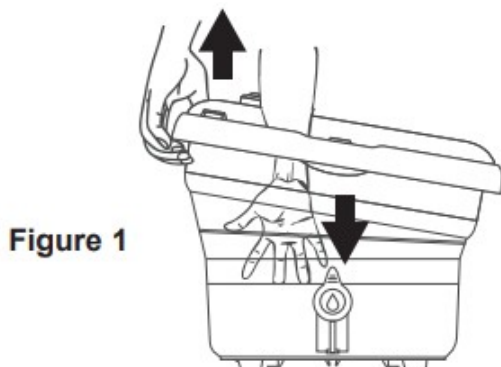
Figure 1

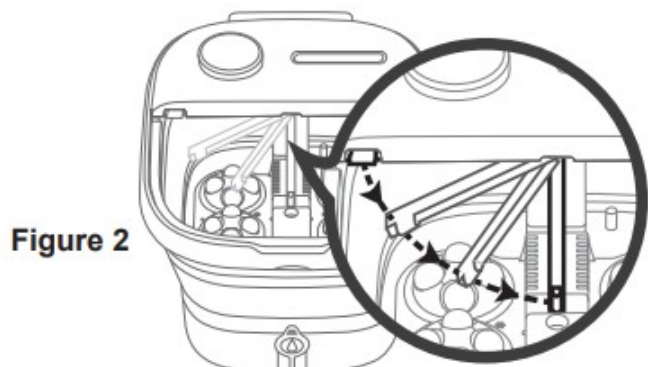
Figure 2

Remove the cover, then keep pressing the machine bottom with one hand, the other lifting one side of the machine body (as shown in Figure 1), then lift the other side in the same way. Press down to install support rod into the bottom

groove(as shown in Figure 2), complete the unfolding.

## FOLD

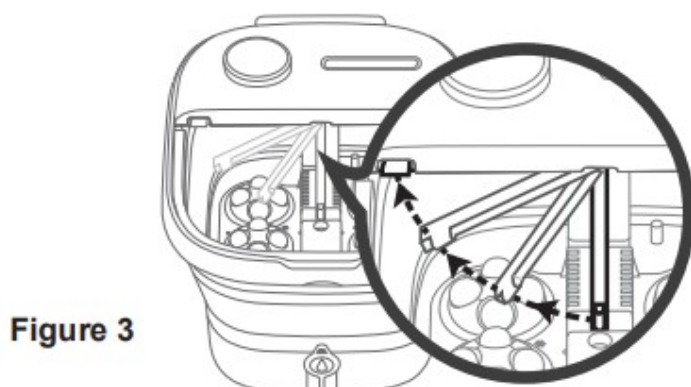
Figure 3

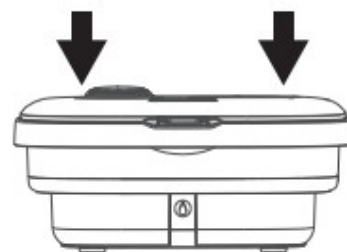
Figure 4

Restore the support rod to its original position (as shown in Figure 3). Simultaneously using both hands to press the machine from the middle position of the top to the state (as shown in Figure 4), complete the folding.

## WARM REMINDER

- Please do not use this machine in the bathroom or other places where the air is humid. If it is inadvertently immersed in water or wet, please do not use it, and it should be handed over to our company or distributor for maintenance and treatment.
- Please use this product in a temperature of 0-40 degrees Celsius.
- Do not store or use the unit near heaters or electric stoves, and keep it away from sunlight.
- This product is for home use only not for commercial use.
- Smells burnt
- Sometimes energized and sometimes de-energized after touching the wire
- Power cord heating
- Other anomalies

## WARNING

- If any of the above phenomena occur, please stop using the product in order to prevent malfunctions and accidents.
- Turn off the power switch, unplug the power plug, and contact the distributor or the company from which you purchased this product.
- Do not disassemble and repair yourself.


## SPECIAL ANNOUNCEMENT

This product has passed strict waterproof test before leaving the factory. There may be a small amount of residual water, which is a normal phenomenon and does not affect product quality. Please feel free to use it!

## WARM REMINDER:

The bottom of this product is designed with a drainage outlet, which can basically empty the residual water in the internal pipe of the machine. However, it is still recommended to wash the interior with clean water after use to prevent the growth of bacteria.

## MASSAGE ROLLER CLEAN

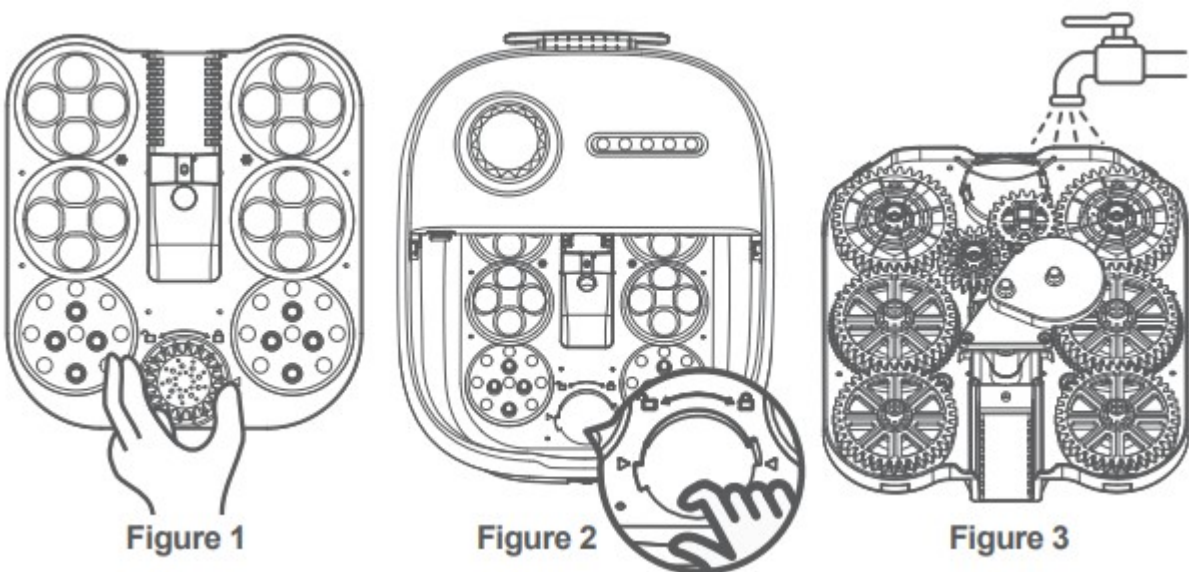First, unscrew the medicine box (Figure 1), then put your finger into the hole shown in Figure 2, and pull the massage ray upward to remove it for easy clearing. For cleaning put the massage tray back in place and press down to complete the installation. If it can not be smoothly put back in place, rotate the massage role by hand and align with the hole to put it back in place.

## SPECIFICATION PARAMETER

|                                 |       |                        |               |
|---------------------------------|-------|------------------------|---------------|
| <b>Model</b>                    | D5    | <b>Model</b>           | FOOT SPA      |
| <b>Rated Voltage</b>            | 220V~ | <b>Max temperature</b> | 48°C          |
| <b>Rated frequency</b>          | 50Hz  | <b>Operation mode</b>  | Touch control |
| <b>Rated power</b>              | 500W  | <b>Display mode</b>    | LED indicator |
| <b>Thermal efficiency grade</b> | >60%  | <b>Heating mode</b>    | PTC heating   |

## COMMON FAULT HANDLING

| Common Fault   | Reason  | Methods of elimination  |
|----------------|---|---|
| No electricity | Power plug unplugged, function switch not turned on | Plug in the power supply in the correct way and turn on the desired switch. If the above methods do not solve the problem, please contact your dealer |
| E1             | Water temperature too high over 53 °C               | When the water temperature drops below 53°C. Please unplug the power cord and reinsert it to restore the machine to normal.                           |
| E2             | Thermostat failure.                                 | Contact seller for replacement and repair   |

## PRECAUTIONS

- If you need to add solid medicine when bathing your feet, please wrap the medicine with gauze and put it in the medicine box to prevent the medicine residue from blocking the filter and the internal pipeline of the machine, causing machine failure.
- It is strictly prohibited to turn on the heating function without water, as it may damage the machine or cause danger!
- It is prohibited to use this machine standing up to avoid damage to the machine and danger!
- Do not operate this machine with wet hands, and do not use this machine near open flames.
- Do not use this unit in wet places such as bathrooms or swimming pools.
- When the ambient temperature is below 0°C, it is strictly forbidden to store water in the product to avoid freezing and damaging the parts.
- If you feel uncomfortable during use (such as dizziness, headache, etc.), please stop using it.
- Do not move the machine during use to avoid fire.
- If liquid leaks from the appliance, the appliance must no longer be used.
- If you find that the machine is malfunctioning or abnormal, you should stop using it immediately and consult or repair it with our company or our company's service centres. Do not disassemble and repair it yourself.
- If the power cord is damaged, it must be replaced by the manufacturer or its maintenance department, or by specialized personnel to avoid danger.

- The machine has a heating function, so it is prohibited for people who are sensitive to heat.
- After use, the switch should be turned off and the power plug should be pulled out of the socket.

### **Disposal of electrical equipment**

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

### **CE declaration**

Malis B. Machoński sp. k. Declares that the Maltec device Water Foot Massager PRO-Hydro500W meets all the requirements of the LVD Directive 2014/35/EU, EC 2014/30/EU and meets the following standards:

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